

# Trends in Amplification

<http://tia.sagepub.com/>

---

## **A Systematic Review of Electric-Acoustic Stimulation : Device Fitting Ranges, Outcomes, and Clinical Fitting Practices**

Paola V. Incerti, Teresa Y. C. Ching and Robert Cowan

*TRENDS AMPLIF* 2013 17: 3

DOI: 10.1177/1084713813480857

The online version of this article can be found at:

<http://tia.sagepub.com/content/17/1/3>

---

Published by:



<http://www.sagepublications.com>

**Additional services and information for *Trends in Amplification* can be found at:**

**Email Alerts:** <http://tia.sagepub.com/cgi/alerts>

**Subscriptions:** <http://tia.sagepub.com/subscriptions>

**Reprints:** <http://www.sagepub.com/journalsReprints.nav>

**Permissions:** <http://www.sagepub.com/journalsPermissions.nav>

>> [Version of Record](#) - Mar 28, 2013

[What is This?](#)

# A Systematic Review of Electric-Acoustic Stimulation: Device Fitting Ranges, Outcomes, and Clinical Fitting Practices

Trends in Amplification  
17(1) 3–26  
© The Author(s) 2013  
Reprints and permissions:  
sagepub.com/journalsPermissions.nav  
DOI: 10.1177/1084713813480857  
tia.sagepub.com  


Paola V. Incerti<sup>1,2,3</sup>, Teresa Y. C. Ching<sup>1,2</sup>, and Robert Cowan<sup>2,3</sup>

## Abstract

Cochlear implant systems that combine electric and acoustic stimulation in the same ear are now commercially available and the number of patients using these devices is steadily increasing. In particular, electric-acoustic stimulation is an option for patients with severe, high frequency sensorineural hearing impairment. There have been a range of approaches to combining electric stimulation and acoustic hearing in the same ear. To develop a better understanding of fitting practices for devices that combine electric and acoustic stimulation, we conducted a systematic review addressing three clinical questions: what is the range of acoustic hearing in the implanted ear that can be effectively preserved for an electric-acoustic fitting?; what benefits are provided by combining acoustic stimulation with electric stimulation?; and what clinical fitting practices have been developed for devices that combine electric and acoustic stimulation? A search of the literature was conducted and 27 articles that met the strict evaluation criteria adopted for the review were identified for detailed analysis. The range of auditory thresholds in the implanted ear that can be successfully used for an electric-acoustic application is quite broad. The effectiveness of combined electric and acoustic stimulation as compared with electric stimulation alone was consistently demonstrated, highlighting the potential value of preservation and utilization of low frequency hearing in the implanted ear. However, clinical procedures for best fitting of electric-acoustic devices were varied. This clearly identified a need for further investigation of fitting procedures aimed at maximizing outcomes for recipients of electric-acoustic devices.

## Keywords

electric and acoustic stimulation (EAS), cochlear implants (CI), hearing aids (HA)

## Introduction

Continuing developments in implant electrode design and improved surgical techniques have resulted in increasing incidence of preservation of residual acoustic hearing in the implanted ear following cochlear implantation (Friedland & Runge-Samuels, 2009; Talbot & Hartley, 2008). In turn, this has led to an expansion of CI candidature to include patients with significantly more residual hearing, and the development of cochlear implant technology for combined electric and acoustic stimulation (EAS) in the implanted ear. In the EAS or Hybrid application, a cochlear implant electrode array is implanted into the cochlea and electrical stimulation is used to convey high frequency information to the user. This is coupled with a hearing aid in the implanted ear, that is used to convey low frequency information to the user via acoustic stimulation.

To date, the focus of EAS research has been on the degree of postoperative hearing preservation outcomes (Arnoldner et al., 2010; Gantz & Turner, 2003; Gstöettner et al., 2006; James et al., 2006; Kiefer et al., 2005; Lenarz et al., 2006;

Mukerjee et al., 2012; Skarzynski, Lorens, Piotrowska, & Skarzynski, 2010; von Ilberg et al., 1999; Woodson, Reiss, Turner, Gfeller, & Gantz, 2009) as well as the perceptual benefits of combined electric-acoustic applications for: speech perception (Büchner et al., 2009; Dorman & Gifford, Dorman, & Brown, 2010; Fraysse et al., 2006; Gantz, Turner, Gfeller, & Lowder, 2005; Helbig & Baumann, 2009; James et al., 2005; Lenarz et al., 2009; Lorens, Polak, Piotrowska, & Skarzynski, 2008; Simpson, McDermott, Dowell, Sucher, & Briggs, 2009; Skarzynski et al., 2012; Turner, Gantz, Karsten, Fowler, & Reiss, 2010); localization (Dunn, Perreau,

<sup>1</sup>National Acoustic Laboratories, Sydney, Australia

<sup>2</sup>The HEARING Cooperative Research Centre, Melbourne, Australia

<sup>3</sup>Department of Audiology, Hearing and Speech Sciences, The University of Melbourne, Australia

## Corresponding Author:

Paola V. Incerti, Maud, National Acoustic Laboratories, Australian Hearing Hub, 16 University Avenue, Macquarie University, NSW 2109 Australia.  
Email: Paola.Incerti@nal.gov.au

Gantz, & Tyler, 2010); perception of music (Brockmeier et al., 2010; Dorman, Gifford, Spahr, & McKarns, 2008; Gfeller, Olszewski, Turner, Gantz, & Oleson, 2006; Gfeller et al., 2007; Gifford, Dorman, & Brown, 2010); and functional performance (Driver & Stark, 2010; Gstoettner et al., 2008; Gstoettner et al., 2011; Helbig et al., 2011). In general, results with EAS have been compared with either the cochlear implant used in isolation, or with the preoperative condition with hearing aids.

A variety of EAS fitting approaches have been reported in the literature. EAS outcome studies report utilising a range of practices for both electric and acoustic stimulation. As the number of recipients increases, research into how to best fit EAS devices to individuals with different degrees of postoperative residual hearing is crucial.

Previously, Talbot & Hartley (2008) have reported a descriptive review of the effectiveness of EAS intervention as compared with conventional cochlear implant use with respect to pitch perception. However, that review was not structured to address three critical issues:

- what is the range of acoustic hearing in the implanted ear that can be effectively used for an electric-acoustic fitting?;
- what benefits are provided by combining acoustic stimulation with electric stimulation?; and
- what clinical fitting practices have been developed for devices that combine electric and acoustic stimulation?

The examination of these issues should lead to the development of guidelines and recommendations for fitting EAS devices.

## Method

### Search Strategy

A systematic review of literature was conducted using methods that followed the guidelines provided by Cox (2005). The initial electronic databases searched included PubMed and MEDLINE. The following key words were entered into the search fields: "Cochlear implant OR cochlear implantation AND acoustic stimulation AND electric stimulation" Additional terms used in subsequent searches: "Cochlear Implantation OR Cochlear Implants AND Electric Stimulation AND Acoustic Stimulation AND Hearing Loss, Sensorineural OR Hearing Loss, High-Frequency. In addition, a textbook by Van de Heyning and Kleine Punte, 2010 was searched by hand to identify references that met the search criteria. Finally, a further electronic search of Web of Science (ISI), SCOPUS - V.4 (Elsevier), CINAHL PLUS (EBSCO), PsycINFO (CSA) databases and the reference lists of the relevant manuscripts were examined for articles that did not appear in prior searches.

### Selection Procedure

The search strategy yielded 168 articles. Following a preliminary examination of all 168 article abstracts, 58 were selected for more comprehensive review. The full text articles were retrieved for the 58 articles and were analyzed with regard to the inclusion and exclusion criteria below.

A final 27 met the criteria and were identified for detailed analysis and included in this review.

### Inclusion Criteria

Studies were selected for detailed analysis if they met the following inclusion criteria:

- publication of the results appeared in a peer-reviewed article or textbook. Excluded repeat publications by the same author/research group using the same subject group and reporting no additional or new evidence.
- randomized controlled trials (RCT), nonrandomized controlled trials, cohort studies, case-control studies;
- studies that implanted "Hybrid or EAS" candidates with preoperative residual low frequency hearing (Pure tone thresholds  $\leq 65$  dB HL for 250 Hz and 500Hz) and severe to profound hearing loss in the high frequencies (Pure tone thresholds  $>70$  dB HL at 1500 Hz);and
- studies that report on preservation or performance data with a minimum of 6 months postoperative device experience only, to ensure the low frequency air-bone gap and cochlear implant MAP/program have stabilized (Helbig et al., 2011; Lenarz et al., 2009).
- articles published between 2000 and 2012

### Exclusion Criteria

Studies were excluded from detailed analysis on the following criteria:

- expert opinion-based evidence and case studies;
- studies that implanted only "traditional" or "conventional" cochlear implant candidates with severe to profound, sensorineural hearing loss;
- studies reporting hearing preservation or performance data with less than 6 months postoperative device experience.

### Rating of Quality Evidence

Two evidence-based practice (EBP) review rating schemes were used to evaluate the research evidence. The first scheme was the Valente et al. (2006) scale design which ranks studies according to the strength of research evidence. This method

was used initially to filter the large number of studies based on: evaluation of the research design used (Level of evidence: 1-6); the quality and relevance of the data (Grade of recommendation: A-D); and the effectiveness (EV: real-world based) and efficacy (EF: laboratory based) of the study. Based on these, an overall strength of recommendation (Strength of evidence: I-III) was made. This process yielded similar strength of evidence for all the studies that met the search criteria. The second scheme was designed by MacDermid (2004). A number of research design principles (subject selection, methodology, endpoints, methods of analysis, results, and conclusions) were used to rate the source of evidence. A total of 24 items are assessed and a total score from (1-48) is calculated. A higher score on the MacDermid scale reflects a higher level of quality of evidence or high internal validity. Whilst a study by Olson and Shinn (2008) showed a strong correlation between the two rating scales used, MacDermid's critical review worksheet "Evaluation of effectiveness Study design" worksheet and "Evaluation Guidelines for Rating the Quality of an Intervention Study" were used together to qualitatively score and rank the evidence from the final 27 studies in more detail. The evidence was also categorized according to whether assessments were carried out under laboratory/ ideal conditions (Efficacious, EF) or in the real world (Effectiveness, EV).

## Results

A total of 27 articles were identified for detailed analysis and inclusion in this review. All articles were published between 2005 and 2012. The review of the studies yielded evidence of a moderate strength (B to D). No meta-analysis or truly randomized controlled trial studies were identified.

## Discussion

### Q1. What is the range of acoustic hearing in the implanted ear that can be effectively used for an electric-acoustic fitting?;

Preservation of residual hearing for the purpose of combining electric and acoustic stimulation has been reported by several research groups using a variety of standard and research electrode arrays and atraumatic surgical techniques (Friedland & Runge-Samuels, 2009; Talbot & Hartley, 2008). To determine the postoperative fitting range in the implanted ear that can be effectively aided, these published studies were classified into three sub-groups: firstly, Type I studies that implanted conventional, perimodiolar electrode arrays using the "Advanced Off Stylet" technique and a cochleostomy approach to achieve a full insertion (300°-430°) of array; secondly, Type II studies that implanted standard, full length or medium length electrode arrays with reduced insertion depth (360°), using either a cochleostomy or round window approach; and

lastly, Type III studies that implanted shorter and thinner electrodes, such as the "Hybrid or EAS arrays," that were specifically designed to increase the likelihood of preserving residual hearing, using either a cochleostomy or round window approach. A summary of the 14 studies that met the search inclusion criteria is shown in Table 1.

The preservation of residual hearing was reported to range from 70% to 100% of participants across the 14 studies with minimum postoperative timeframes of 6 to 12 months. In these studies the term "Hearing preservation" encompassed both complete and partial preservation. As defined in the studies complete preservation referred to the maintenance of hearing thresholds within 10 dB HL of all preoperative thresholds, whereas partial preservation generally referred to a change in thresholds of greater than 10 dB to 15 dB HL of any preoperative threshold level. The degree of postoperative hearing loss and audiometric configuration varied both across individuals within studies, and across studies (see Table 1: Key Findings & Range).

No systematic reporting procedure exists to allow for direct comparisons of postoperative hearing results across various electrode arrays or surgical procedures aimed at improving hearing preservation. In order to compare changes in low frequency acoustic hearing across Type I, II & III studies, the preoperative and postoperative medians at 125 Hz, 250 Hz & 500 Hz were calculated using the raw data from each participant's audiogram, for six studies out of 14 that provided individual threshold data across frequencies (Gstoettner et al., 2006, 2009; James et al., 2005; Podskarbi-Fayette, Pilka, & Skarzynski, 2010; Skarzynski et al., 2012; Simpson et al., 2009). Thresholds at  $\geq 750$  Hz were excluded because reporting of thresholds at the octave and intermediate frequencies was highly variable between studies and the preoperative thresholds were already in the severe to profound range. Figure 1 shows the median change in low frequency hearing at 125, 250, and 500 Hz. Results from each study are represented by a data point  $\pm$ SD for the various electrode array studies.

Figure 1 shows that the greatest shift in acoustic hearing thresholds occurred at 500 Hz across all studies. Median changes in postoperative hearing thresholds ranged from 5 dB to 25 dB HL at 125 Hz, 10 dB to 30 dB HL at 250 Hz, and 28 to 43 dB HL at 500 Hz. When examining threshold shifts it is important to also consider a shift in the context of the impact on individual hearing thresholds and whether they are still aidable. Threshold shifts may move postoperative hearing thresholds into a mild to moderate range that may still benefit from amplification or may also result in final thresholds that are in the severe to profound range that may not.

Table 2 summarizes the hearing preservation data for Type I, II & III studies based on six of the 14 studies that reported individual threshold data across frequencies.

Average changes to thresholds at and below 500 Hz (3FA) ranged between approximately 24 dB to 37 dB HL and median changes ranged from 18 dB HL to 34 dB HL. In

**Table I.** A Summary of 14 Articles That Provide Information on Postoperative Hearing Preservation Outcomes That Met the Systematic Review Inclusion Criteria.

Type I. Perimodiolar electrode array ( "Advance off styllet" technique) studies					
Study	Study design	Implant/electrode/surgical approach	Key findings	Level, grade strength, EV/EF	MD score
James et al., 2006	I group: N = 10 adults	Implant: Nucleus 24 Contour Advance. Electrode array: 19mm TL, 22 channel Cochleostomy (<1.5mm) Anterior & inferior to RW. Insertion depth to 17-19mm or (300°-430°).	70% preservation of hearing to varying degrees for participants at minimum of 6 months. Range: Individual thresholds in implanted ear from normal to profound at 125-500Hz & profound at ≥750Hz.	3B II EF	25/48
Simpson, McDermott, Dowell, Sucher, & Briggs, 2009	I group: N = 5 adults	Implant: Nucleus Freedom Contour Advance Electrode array: 22 channel, 9mm TL. Cochleostomy: Inferior to RW. "AOS" technique. Insertion depth to 17-19mm or (300°-430°).	100% preservation to varying degrees for trial participants after a minimum of 6 months. Range: Individuals thresholds in implanted ear from mild to profound HL range at 125-500Hz & profound at 750Hz.	3B II EF	30/48
Type II. Standard length (Reduced/shallow insertion depth) and/or medium length straight electrode array studies					
Study	Study design	Implant/electrode/surgical approach	Postoperative hearing preservation	Level, grade strength, EV/EF	MD score
Kiefer et al., 2005	I group: N = 13 adults	Implant: MED-EL Combi 40+ Electrode array: 12 channel, 1.Standard/flex 31.5mm TL (n = 6) 2.Medium, 26.4mm TL (n = 7). Cochleostomy. Insertion depth to 18-24mm or 360°	84.6 % preserved hearing either fully or partially reported at minimum of 12 months. Range: No individual data reported at 12 months.	3C II EF	24/48
Gstoettner et al., 2006	I group: N = 23 adults	Implant: MED-EL Combi 40+ Electrode array: 12 channel 1.Standard/flex array 31.5mm TL 2.Medium array, 26.4mm TL. Cochleostomy. CT scans to predict insertion depth to achieve 360° (18-24mm). Reduced insertion for longer arrays.	70% preservation to varying degrees for trial participants after a minimum of 6 months (range 6 -70 months). Range: Individuals thresholds in implanted ear from normal to profound range at 125-500Hz & severe to profound at 750Hz.	4C II EF	20/48

(continued)

**Table I. (continued)**

Type II. Standard length (Reduced/shallow insertion depth) and/or medium length straight electrode array studies					
Study	Study design	Implant/electrode/surgical approach	Postoperative hearing preservation	Level, grade strength, EV/EF	MD score
Skarzynski, Lorens, Piotrowska, & Anderson, 2007	1 group: N = 9 children Inclusion criteria: "partial deafness"	Implant: MED-EL Combi 40+ or Pulsar. Electrode array: 12 channel 1. Standard/flex 31.5mm TL (n = 1) 2. Medium, 26.4mm TL (n = 8) RW approach. Reduced insertion depth (20mm) for longer arrays	100% preservation either partially or fully for trial participants after a minimum of 6 months (range 6-12 months). Range: Individuals thresholds in implanted ear from normal to severe HL range at 125-500Hz & moderate to profound at 750Hz .	4C II EF	23/48
Gstoettner et al., 2008	1 group: N = 18 adults	Implant: MED-EL Combi 40+ Electrode array: 12 channel Medium array, 26.4mm TL Cochleostomy: Limited electrode insertion of 18mm-22mm.	83.2 % preserved hearing either fully or partially for trial participants after a minimum of 12 months. Range: No individual data reported at 12 months.	4C II EF & EV	24/48
Skarzynski, Lorens, Piotrowska, & Podskarbi-Fayette, 2009	2 groups: Total =28 N = 19 adults N = 9 children	Implant: MED-EL Combi 40+ or Pulsar. Electrode array: 12 channel 1. Standard/flex, 31.5mm TL (n = 18) 2. Medium, 26.4mm TL,(n = 10). RW approach. Reduced insertion depth for longer arrays (20mm).	84% preserved hearing for trial participants either partially or fully after a minimum of 12 months (range 1-4 years postoperatively). Range: No individual data. Means reported.	3C II EF	26/48
Type III. Electroacoustic electrode array studies					
Study	Study design	Implant/electrode/surgical approach	Postoperative hearing preservation	Level, grade strength/EF	MD score
Gantz, Turner, Gfeller, & Lowder, 2005	1 group N = 21 adults	Implant: Nucleus Hybrid 10. Electrode array: 6 channel. Short array 10 mm (T L) Cochleostomy (0.5mm)	95.8 % preserved hearing either partially or fully after minimum of 6 months (range 6 months - 36 months). Range: No individual data. Means reported.	3C II EF	24/48
Gstoettner et al., 2009	1 group N = 9 total	Implant: MED-EL Pulsar Electrode array: 12 channel Flex EAS 6.4mm T L RW approach (n = 7) Cochleostomy (n = 2). Insertion depths to achieve 360° (18-22mm). Calculated to enter the 1000 Hz region.	100 % preserved hearing either partially or fully after minimum of 6 months (range 6 - 15.65 months). Range: Individuals thresholds in implanted ear from normal to profound HL range at 125-500Hz & severe to profound at 750Hz.	4C II EF	22/48

(continued)

Table 1. (continued)

Type III. Electroacoustic electrode array studies		Study design	Implant/electrode/surgical approach	Postoperative hearing preservation	Level, grade strength/EF	MD score
Lenarz et al., 2009	2 groups N = 32 adults 1. Hybrid candidates 2. Conventional CI candidates	Implant: Nucleus Hybrid—L24 Electrode array: 22 channel L array, 18mm TL RW approach & and full insertion (depth of 16mm or 270°)	100 % preserved hearing. * * PTA results reported at * 3 months (n = 32/32) 6 months (n = 22/32) & 12 months (n = 16/32) postoperatively. Average increase in air-bone gap was 4.4 dB at 1 month. At 6 months there was an average increase of 1.2 dB compared to preoperative value.	3C II EF	27/48	
Danka, Pillsbury, Adunka, & Buchman, 2010	2 groups: N = 20 adults EAS array (n = 10) & Standard CI array (n = 10) groups matched for age & word scores.	Implant: MED-EL Combi 40+ & Pulsar CI100. Electrode array: 12 channel 1. Flex EAS, 26.4 mm TL (n = 10) 2. Standard straight 31mm TL full insertion (n = 10) EAS group: RW approach (n = 5) & Cochleostomy (n = 5)	Postoperative hearing preservation: minimum of 6 months device use. EAS array: 90 % preserved hearing Standard array: 0% preserved hearing Range: No individual data. Means reported	2B II EF	32/48	
Woodson et al., 2009	1 group N = 81 total Phase I & II: <60 dB HL below 500Hz	Implant: Nucleus Hybrid 10 Electrode array: Short array, 6 channel, 10 mm (TL) Cochleostomy (0.5mm) anterior & inferior to RW.	90 % preserved hearing either partially or fully after minimum of 12 months (ranging from 12 months to 5 years). Range: Individuals thresholds in implanted ear from normal to profound HL range at 125-750 Hz.	3B II EF	29/48	
Helbig et al., 2011	1 group N = 18 adults	Implant: MED-EL Pulsar CI100 Electrode array: Flex EAS, 12 channel, 26.4 mm TL RW insertion or Cochleostomy	100 % preserved hearing either partially or fully after minimum of 15 months postoperatively.	3B II EF & EV	28/48	
Skarzynski et al., 2012	3 group N = 23 adults 1. EC group ≤ 50 dB HL at 500Hz (n = 7) 2. EAS group Between 50 and 80dB HL (n = 7) 3. ES group ≥ 80 dB HL.	Implant: Nucleus Straight Research Array Electrode array: 22 channel, 25 mm TL RW insertion	Hybrid candidates PTA results reported at * 13 months postoperatively (EC group n = 7/7 & EAS group 5/7 EC). 100 % preserved hearing either partially or fully. Slight improvement between 1 month and 4 months postoperatively due to clearing of middle ear fluid.	3C II EF	27/48	

## Abbreviations:

AOS= Advanced Off Stylet; CI= Cochlear implant; EAS=Electric-acoustic Stimulation; EC=Electrical complement; EF=Effacious or laboratory based measures; ES=Electrical stimulation; EV= Effectiveness or real world measures; IE= Implanted ear; HA= Hearing aids; MD score= MacDermid score; NIE= No implanted ear; NH: Normal hearing; PDCI; Partial deafness cochlear implantation; PTA: pure tone audiometry; PTT: pure tone thresholds; RW= round window; SN= Sensorineural hearing loss; SNR = signal to noise ratio; TL=Total length.

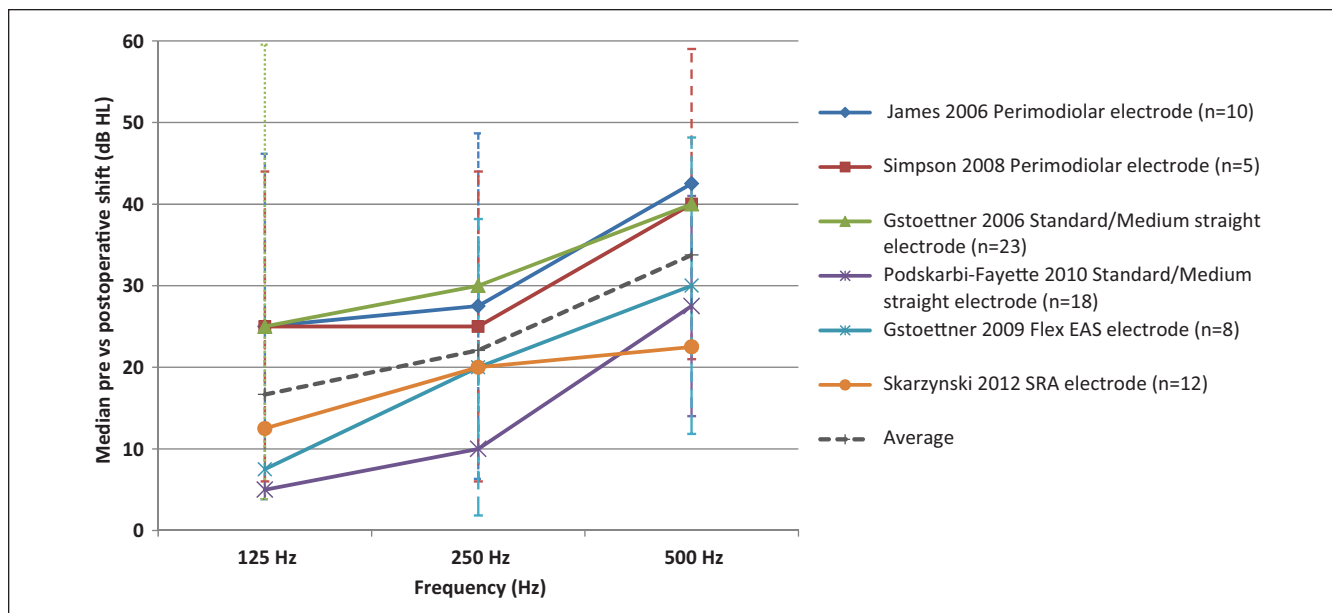


Figure 1. Median change in low frequency thresholds dB HL ± SD for various electrode array studies.

Table 2. Hearing Preservation Data for Type I, II & III Electrode Array Studies.

Array	Total number participants	Average threshold shift in dB HL 3FA (125,250 & 500 Hz)	Median threshold shift in dB HL & 25th, 75th percentile (125, 250 & 500 Hz.)			Postoperative Pure tone thresholds minimum –maximum dB HL			Complete preservation <10dB HL at 3 frequencies	Partial preservation ≥15dB HL at 1 frequency or ≥10dB HL at 2 adjacent frequencies	Total loss No measurable hearing thresholds
			Median	25th	75th	125 Hz	250 Hz	500 Hz			
Type I James et al. <sup>1</sup> (2006) Simpson, McDermott, Dowell, Sucher, & Briggs <sup>2</sup> (2009)	17	37 <sup>1,2</sup>	34 <sup>1,2</sup>	25 <sup>1,2</sup>	45 <sup>1,2</sup>	10-90	25-110	65-125	0 (0%)	14 (82%)	3 (18%)
Type II Gstoettner et al. <sup>1</sup> (2006) Podskarbi-Fayette, Pilka, & Skarzynski <sup>2</sup> (2010)	41	26 <sup>1,2</sup>	22 <sup>1,2</sup>	17 <sup>1,2</sup>	38 <sup>1,2</sup>	20-90	25-110	45-125	13 (32%)	21 (51%)	7 (17%)
Type III Gstoettner <sup>1</sup> (2009) Skarzynski <sup>2</sup> (2012) Helbig et al. <sup>3</sup> (2011)	26	24 <sup>1,3</sup>	18 <sup>1,2</sup>	18 <sup>1,2</sup>	23 <sup>1,2</sup>	15-50 <sup>1</sup>	30-90 <sup>1</sup>	70-120 <sup>1</sup>	2 <sup>1</sup> (22%)	7 <sup>1</sup> (78%)	0 <sup>1</sup> (0%)

Note: <sup>1,2,3</sup>Super script refers to the relevant reference in the Type row the calculation was based on.

general, better preservation outcomes were achieved in Type III studies that implanted electrode arrays specifically designed to preserve residual hearing for the purpose of electric and acoustic stimulation compared to Type I and II studies that implanted standard electrode arrays either partially or fully inserted. Table 2 shows a lower average threshold shift (3FA) postoperatively for Type III studies as compared to either Type I or Type II. In addition, a higher percentage of participants showed complete preservation (defined in this paper as <10dB HL at any 3 frequencies) or partial preservation (defined as ≥15dB HL at 1 frequency or ≥10dB HL at any 2 adjacent frequencies) of residual hearing for Type III

studies as compared to either Type I or Type II . Finally, Type I and Type II studies reported total hearing loss percentages of 17% and 18 % respectively, whereas, Type III array studies report no cases of a participant losing all hearing as a result of surgery. There is a wide range of resulting auditory thresholds in the low frequencies following implantation (see Table 2; minimum and maximum pure tone thresholds). Low frequency thresholds in the implanted ear (125, 250, and 500 Hz) ranged from normal to profound. However, the majority of subjects showed thresholds in the severe to profound range for frequencies of 750 Hz and above. The data are too limited at this time to provide a definitive recommendation regarding



the risk of this intervention with regard to the various electrode arrays.

It is important to note that the majority of studies that met the inclusion criteria had relatively small numbers of participants and reported results with minimum postoperative timeframes of 6 to 12 months. Whilst a few studies reported outcomes for longer timeframes, unfortunately, these results were combined with results from much shorter time periods. Examination of the study evidence revealed a high degree of variability in assessment time frames and reporting procedures, therefore, the evidence is insufficient to provide any recommendation regarding the longer term changes in postoperative hearing following this intervention.

In summary, studies investigating preservation of residual acoustic hearing have shown that preservation of low frequency hearing was achieved in the majority of patients, particularly with the use of electric-acoustic arrays. However, there is on average, a reduction in the low frequency acoustic thresholds postoperatively. On an individual level, there is a wide range of resulting auditory thresholds and audiometric configurations, with the majority available for EAS application following implantation. Given this, there is a need for a range of amplification options to accommodate degrees of hearing loss from mild to profound in the implanted ear. The review also emphasizes the need for additional research studies utilizing longer timeframes with minimum of 2 to 5 years postoperatively and larger studies as postoperative hearing preservation is not always stable and changes in postoperative hearing can continue to occur. It is recommended that studies do not combine hearing preservation data from various postoperative timeframes or exclude any enrolled participants. Standard reporting of thresholds at octave and intermediate frequencies, median and mean postoperative threshold changes at all frequencies and preoperative and postoperative audiometric data for both ears for all enrolled participants would enable meta-analysis methods to be used to examine preservation results across studies, across surgical technique and electrode arrays.

A systematic and consistent reporting procedure for all hearing preservation results would also help clinicians in the preoperative counseling stage to provide a clearer understanding of risks of this intervention and assist patients considering this option in making a fully informed choice. An international consensus statement on recommended reporting practices for postoperative hearing preservation and performance data for EAS devices would be clinically useful and relevant.

## **Q2. What benefits are provided by combining acoustic stimulation with electric stimulation?**

Following review, 18 studies were selected for detailed analysis to address the second question concerning benefit of combining acoustic and electric stimulation. Findings are summarized in Table 3.

## **Preoperative Versus Postoperative Speech Perception**

To evaluate the speech perception benefits resulting from combining acoustic amplification with electric stimulation, preoperative speech perception scores obtained using hearing aids (HAs) were compared to postoperative scores in the implanted ear using electric stimulation alone (CI alone) and combining acoustic and electric stimulation (CI +HAI). Results consistently show that speech perception, on average, was significantly higher postoperatively with electric stimulation alone (CI) when compared to best aided preoperative condition, and for combined acoustic and electric stimulation in the implanted ear (CI+HAI) as compared to scores for electric stimulation alone (CI). These findings are shown in Figure 2 for the various electrode arrays studies. In all studies, speech was presented from a single frontal loudspeaker (0° azimuth) at a fixed distance from the participant. Word perception scores from listening in quiet are shown, with averaged results from each study calculated and represented by a data point.

Figure 2 shows speech performance in the preoperative best-aided condition (HAs) compared with postoperative monaural listening scores obtained for the implanted ear using electric stimulation alone (CI only) represented as the closed symbols. Any increase in measured performance is assumed to represent primarily benefit as a direct result from implantation. The graph shows that all studies, regardless of array used, reported significantly better mean word perception scores for postoperative electric stimulation alone (CI only) as compared to preoperatively (HAs) in quiet. The increase in word scores ranged from 30% to 39% across the various electrode studies. It could be argued, that testing this condition provides evidence to support this intervention for individuals with severe high frequency sensorineural hearing losses, even when there is complete hearing loss associated with the implantation.

Figure 2 also compares speech performance in the preoperative best aided condition (HAs) as compared with the unilateral listening condition in the implanted ear using combined acoustic and electric stimulation (CI+HAI). Data was taken from the same studies and represented as open symbols. Any increase in measured post implantation performance represents primarily benefit from use of preserved low frequency residual hearing through combining acoustic and electric stimulation following cochlear implantation. On average, word scores in quiet were significantly better postoperatively in the implanted ear combining acoustic and electric stimulation (CI+HAI) as compared to preoperatively (HAs). Postoperative EAS benefit ranged from 43.8% to 49% in word scores. The post implantation speech results show that there is benefit combining acoustic and electric stimulation.

Across all studies eligible for inclusion, results consistently show that speech perception, on average, was enhanced postoperatively with electric stimulation alone compared to

**Table 3.** A Summary of 18 Articles That Provide Information on Postoperative Outcomes That Met the Systematic Review Inclusion Criteria.

1. Perimodiolar electrode array ( "Advance off styllet" technique) studies					
Study	Participants	Outcome measures	Key findings	Level, Grade Strength, EV/EF	MD Score
James et al., 2006	I group: N = 10 adults Inclusion criteria: Word scores $\geq 10\%$ (IE). Symmetrical.	Words at 65 dB SPL in quiet; German, Spanish & French Sentences at 70 dB SPL in multitalker babble (+5dB SNR) Single loudspeaker at 0° azimuth.	CI alone scores significantly better than preoperative scores with HAs, in quiet & noise. (N = 7/10) Mean CI + Ipsilateral HA scores better than CI alone scores in quiet & noise (not statistically significant). Ceiling effects in quiet (N = 2) and in noise (N = 1). Results reported for subjects with HTLs $\leq 80$ dB HL at 125 & 250 Hz in IE.	3B II EF	25/48
Simpson, McDermott, Dowell, Sucher, & Briggs, 2009	I group: N = 5 adults Inclusion criteria: Word scores at 60-65 dBA between 0-30% in both ears	Words at 60-65 dB SPL in quiet; English Sentences at 65 dB SPL in adaptive noise (SNR50). English Single loudspeaker at 0° azimuth.	CI alone word scores significantly better than preoperative scores in quiet. (N = 3/5 bimodal & N = 2/5 CI+ bilateral HAs) CI +HA/s scores significantly better than CI alone scores in quiet & noise. No significant difference between CIHAi, CIHAc, CIHAs condition (N = 2)	3B II EF	30/48
2. Standard length (Reduced/shallow insertion depth) and/or medium length straight electrode array studies					
Study	Participants	Outcome measures	Postoperative findings	Level, grade strength, EV/EF	MD score
Kiefer et al., 2005	I group: N = 13 adults Inclusion criteria: Word scores at 70 dB $\leq 40\%$ in both ears.	Words at 70 dB SPL in quiet Sentences in quiet and noise (+10dB SNR) German Single loudspeaker at 0° azimuth.	CI alone word scores significantly better than preoperative scores in quiet & noise. CI +HA/s sentences scores significantly better than CI alone scores in quiet & noise. Results reported for participants who retained hearing (N = 11/13).	3C II EF	24/48
Gstoettner et al., 2006	I group: N = 23 adults Inclusion criteria Word scores at 70 dB $\leq 40\%$ in best aided condition	Words at 65 dB SPL; Sentences at 65 dB SPL in adaptive noise (SNR50). German Single loudspeaker at 0° azimuth.	Mean CI alone word scores were better than preoperative scores in quiet & noise. No statistical test performed on data. Mean CI + Ipsilateral HA sentence scores were better than CI alone scores in quiet. No statistics. Ceiling effect for words in quiet Results reported for participants with complete hg preservation (N = 8/23).	4C II EF	20/48

(continued)

Table 3. (continued)

Study	Participants	Outcome measures	Postoperative findings	Level, grade strength, EV/EF	MD score
2. Standard length (Reduced/shallow insertion depth) and/or medium length straight electrode array studies					
Lorens et al, 2007	3 groups: Partial deafness adults ( $n = 11$ ) Traditional CI adult recipients ( $n = 22$ ) and NH adults ( $n = 20$ ).	Words in quiet and noise (+10 dB SNR & +0 dB SNR). Polish Sentences in noise (+10dB SNR) Single loudspeaker at 0° azimuth.	EAS group performed significantly better than traditional CI recipients for word & sentences. EAS group: CI + bilateral HAs word scores were significantly better than CI alone, HA alone & CI + Ipsilateral HA scores in noise. CI + Ipsilateral HA word scores were significantly better than CI alone, HA alone in noise. Ceiling effect for words in quiet & sentences in noise.	3B II EF	27/48
Skarzynski, Lorens, Piotrowska, & Anderson, 2007	I group: $N = 9$ children Inclusion criteria: "partial deafness" or significant residual hearing in low frequencies.	Words in quiet and noise (+10 dB SNR) Polish Single loudspeaker at 0° azimuth.	Postoperative word scores* significantly better than preoperative scores with HAs in quiet. ( $n = 4/9$ ). * Different devices used: 4/9 used residual hearing with CI. 4/9 used a DUET combined EAS stimulation device. 1/9 uses CI only.	4C II EF	23/48
Helbig, Baumann, Helbig, von Malsen-Waldkirch, & Gstoettner, 2008	I group: $N = 9$ adults Inclusion criteria: EAS recipients.	Words at 70 dB SPL in quiet. Sentences in quiet and in noise (+10dB, +5dB & 0 dB SNR). German. Single loudspeaker at 0° azimuth. Abbreviated Profile of Hearing Aid Benefit (APHAB) questionnaire	DUET (combined electric & acoustic stimulation device) sentences scores were significantly better than CI alone scores in noise (+10 dB & +5dB SNR). APHAB questionnaire No significant difference reported with device upgrade from CI + Ipsilateral ITE to DUET Ceiling effects for sentences in quiet & floor effects for sentences in noise (0 dB SNR).	3B II EF & EV	24/48
Gstoettner et al., 2008	I group: $N = 18$ adults Inclusion criteria: Word scores $\leq 45\%$ in best aided condition.	Words at 70 dB SPL in quiet; German, Spanish & English Sentences in quiet and noise (+10 dB SNR) Single loudspeaker at 0° azimuth. Abbreviated Profile of Hearing Aid Benefit (APHAB) questionnaire.	CI + Ipsilateral HA word & sentence scores were significantly better than preoperative scores with HAs in quiet & noise. APHAB questionnaire significant postoperative improvement in subjective benefit compared to preoperative HA condition. Speech testing results reported for 6 out of 18 participants who retained hearing and choose to wear amplification.	4C II EF & EV	24/48
Skarzynski, Lorens, Piotrowska, & Podskarbi-Fayette, 2009	$N = 28$ total 2 group $N = 19$ adults $N = 9$ children Inclusion criteria: "partial deafness"	Words at 60 dB SPL in quiet and noise (+10 dB SNR) Polish Single loudspeaker at 0° azimuth.	Word scores in best aided condition* were significantly better than preoperative scores with HAs in quiet & noise. *Best aided condition: DUET + contralateral acoustic hearing or CI + natural bilateral acoustic hearing.	3B II EF	26/48

(continued)

**Table 3. (continued)**

2. Standard length (Reduced/shallow insertion depth) and/or medium length straight electrode array studies				
Study	Participants	Outcome measures	Postoperative findings	Level, grade strength, EV/EF MD score
Podskarbi-Fayette, Pilka, & Skarzynski, 2010	N = 18 total N = 11 adults N = 7 children Inclusion criteria: Residual hearing at levels in the low frequency range to 500Hz-1KHz	Monosyllabic words at 60 dB SPL in quiet and noise (+10 dB SNR). Polish Single loudspeaker at 0° azimuth	Word scores in best aided condition* were better than preoperative scores with HAs in quiet & noise. No statistical test performed on data. (N = 11/18) *Best aided condition: DUET + contralateral acoustic hearing (n = 7) or CI + natural bilateral acoustic hearing (n = 10) or CI alone (n = 1).	4C II EF 23/48
3. Straight electroacoustic electrode array studies				
Study	Participants	Outcome measures	Key findings	Level, grade strength, EV/EF MD score
Gantz, Turner, Gfeller, & Lowder, 2005	I group N = 21 adults Inclusion criteria: Word scores at 70dB SPL between 10-60% (IE) & 80% (NIE)	Words at 70 dB SPL in quiet. Sentences in quiet and noise (+10 dB SNR). English Single loudspeaker at 0° azimuth	CI alone word scores were better than preoperative HA scores in quiet (N = 11/21). Mean CI +HAs SRT scores in noise were better than retrospectively matched group of traditional CI users. No statistics.	3C II EF 24/48
Buchener et al., 2009	I group N = 22 adults Inclusion criteria Hybrid candidates ≤ 60dB up to 500Hz	Sentences at 65 dB SPL in (+10 dB SNR) and adaptive noise (SNR50) Oral articulation evaluation Single loudspeaker at 0° azimuth	CI + Ipsilateral HA scores were significantly better than CI or HA alone in noise. Information up to 300Hz significantly improved speech perception when combined with CI signal compared to CI alone. Significant correlation between preoperative articulation abilities & CI alone SRT scores.	3B II EF 24/48
Gstoettner et al, 2009	I group N = 9 adults Inclusion criteria Word scores at 65 dB SPL in quiet ≤ 40% best aided condition	Words at 65 dB SPL in quiet. Sentences in quiet and noise (+10 dB SNR). German Single loudspeaker at 0° azimuth	CI + Ipsilateral HA scores were higher than CI alone or HA alone scores in quiet and noise. (N = 7/9) No statistics.	4C II EF 22/48
Lenarz et al, 2009	2 groups N = 32 adults 1. Hybrid candidates ≤ 60dB up to 500Hz (n = 24) 2. Conventional CI candidates (n = 8)	Words at 65 dB SPL in quiet. Sentences at 65 dB SPL in adaptive noise (SNR50). German Single loudspeaker at 0° azimuth	CI alone scores were significantly better than preoperative scores in quiet & noise. (N = 22/32) CI + Ipsilateral HA sentence scores were significantly better than CI alone & HA alone in noise. Adults with hearing loss >30 years (n = 6) showed limited benefit in CI alone condition. Hearing loss of <30 years showed large EAS benefit. Substantial inter-individual differences in scores.	3C II EF 28/48


(continued)

Table 3. (continued)

Study	Participants	Outcome measures	Key findings	Level, grade strength, EV/EF	MD score
3. Straight electroacoustic electrode array studies					
Adunka, Pillsbury, Adunka, & Buchman, 2010	2 groups N = 20 total Inclusion criteria Word scores at 65 dB SPL in quiet ≤ 50% in best aided condition	Words in quiet. English Single loudspeaker at 0° azimuth 1. EAS array recipients (n = 10) 2. Standard CI array recipients (n = 10) (Groups matched for age & preoperative word score)	CI alone scores were better than preoperative scores in quiet. (N = 19/20) EAS electrode recipient group: CI + Ipsilateral HA word scores were significantly better than CI alone & HA alone in quiet. The EAS electrode group word scores were significantly better than the conventional CI group in quiet.	2B II EF	32/48
Dunn, Perreau, Gantz, & Tyler, 2010	1 group N = 11 total Inclusion criteria Hybrid recipients.	The Recognition with Multiple Jammers speech-perception-in-noise test. English Everyday sounds localization test. Eight speaker array.	Bilateral acoustic amplification with electric stimulation in implanted ear significantly better for sentence recognition compared to all other conditions (CI alone, HA alone, Bimodal) in spatially separated noise (N = 9/11). Significant difference in localization for the binaural HA conditions compared to monaural hearing aid conditions.	3B II EF	28/48
Gifford, Dorman, & Brown, 2010	4 groups: N = 64 adults 1. EAS array (n = 5) 2. Unilateral CI array (n = 25) 3. Bilateral CI array (n = 10) 4. Bimodal (n = 24)	R_SPACE 8 loudspeakers Adaptive speech test using spatially separated noise. English Frequency resolution at 500Hz by derived auditory filter shapes Words at 65 dB SPL in quiet. Sentences in quiet and noise (+10 dB SNR) German & Dutch Single loudspeaker at 0° azimuth	EAS group (n = 5) combined SRT scores significantly better than Unilateral CI group (n = 25), Bilateral CI group (n = 10) & bimodal (CI + contralateral HA) group (n = 24). Normal hearers have significant better frequency selectivity than preoperative EAS participants (n = 5). CI alone scores were significantly better than preoperative scores in quiet & noise. (N = 17/18)	3B II EF	27/48
Helbig et al., 2011	1 group N = 18 adults Inclusion criteria Word scores at 65 dB SPL in quiet ≤ 50% in best aided condition	Words at 65 dB SPL in quiet. Sentences in quiet and noise (+10 dB SNR) German & Dutch Single loudspeaker at 0° azimuth	CI alone scores were significantly better than preoperative scores in quiet & noise. (N = 17/18) CI + Ipsilateral HA speech scores were significantly higher than CI alone or HA alone scores in quiet & noise.	3B II EF & EV	28/48
Skarzynski et al., 2012	3 group N = 23 adults Hybrid candidates (EC & EAS) and Conventional CI candidate (ES).	Abbreviated Profile of Hearing Aid Benefit (APHAB) questionnaire Words at 60 dB SPL in quiet and noise (+10 dB SNR) Polish Single loudspeaker at 0° azimuth.	APHAB questionnaire significant postoperative improvement in subjective benefit is reported compared to preoperative HA condition. Postoperative word scores were better than preoperative scores with HAs in quiet & noise for "Hybrid" candidate group (EC group n = 7/7 & EAS group 5/7 reported at 13 months). *Postoperative aided device condition not reported.	3C II EF	27/48

Abbreviations:

APHAB=Abbreviated Profile of Hearing Aid Benefit questionnaire; BTE= Behind the ear; CI: Cochlear implant; dB = decibel; EAS= Electric-acoustic Stimulation; EC=Electrical complement; EF=Efficacious or laboratory based measures; ES=Electrical stimulation; EV= Effectiveness or real world measures; HA = hearing aids; IE= Implanted ear; ITE=In the ear; HA= Hearing aids; MD score= MacDermid score; NIE = nonimplanted ear; NH= Normal hearing; PDCI=Partial deafness cochlear implantation; PTA= Pure tone Audiometry; PTT=Pure tone thresholds; SN= Sensorineural hearing loss; SNR = signal to noise ratio.

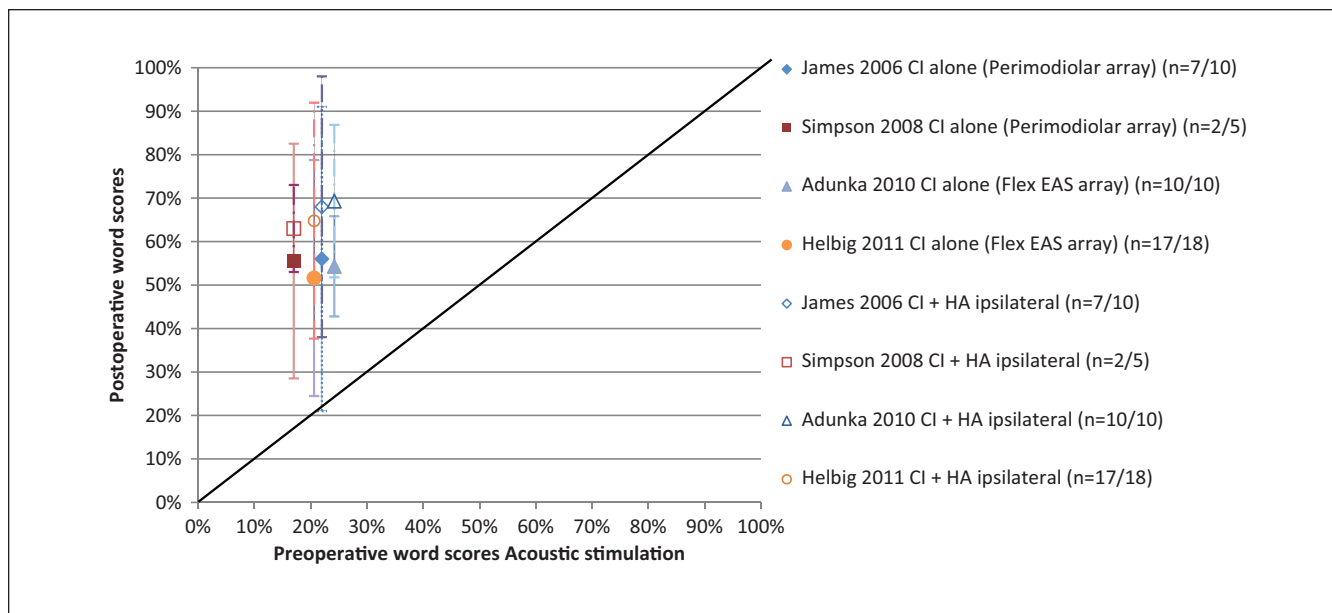

**Figure 2.** Preoperative versus postoperative word scores in quiet (Mean ± SD) in implanted ear for various array studies. Closed symbols represent electric stimulation alone. Open symbols represent combined acoustic & electric stimulation.

preoperative condition. Greater postoperative advantage was evident when combining acoustic and electric stimulation in the implanted ear than for use of electric stimulation alone. It should also be noted that considerable individual variability was reported across these studies, with the majority of individuals showing benefit, whereas a small number showed no difference in preoperative versus postoperative performance. Age and duration of hearing loss accounted for significant variance and had a negative impact on speech performance (Gantz et al., 2009; Lenarz et al., 2009). A significant correlation was also found between preoperative articulation abilities and postoperative CI alone Speech Reception Threshold (SRT) sentence scores (Buchener et al., 2009).

### Acoustic and Electric Stimulation Speech Perception Benefit in the Implanted Ear

A further question is the degree of benefit that can be gained from combining acoustic and electric stimulation when residual hearing is preserved following implantation? To investigate this, several studies have compared speech perception in the implanted ear using electric stimulation only (CI) with speech perception using combined acoustic and electric stimulation (CI+HAi) in quiet (Figure 3a) and in noise (Figure 3b). Speech perception scores for testing in quiet and noise are shown, with averaged results from each study calculated and represented by a data point.


As shown, speech perception scores for the implanted ear were consistently higher for the combined acoustic and electric stimulation as compared to electric stimulation alone, particularly for speech perception tests in noise.

Figure 3a shows that in quiet, benefit ranged from 6% to 15% for word perception scores. Five out of the nine electrode array studies reported statistically significant benefit on average. Figure 3b shows that in noise (+10dB SNR) this benefit ranged from 4% to 29% for sentence perception scores. Five out of six studies reported statistically significant benefit, on average. The one study that did not reach statistical significance had reached ceiling affect.


Across all studies, speech perception scores were higher postoperatively in the implanted ear when using combining acoustic and electric stimulation (CI+HAi) as compared to electric stimulation alone (CI), on average, in both quiet and noise. The addition of acoustic stimulation potentially provides low frequency information to enhance the perception of F0 voicing cues (timing and presence), envelope periodicity cues, low-frequency segmental phonetic cues (nasal resonant, F1 onset and shape of transitions) for improved speech perception (Brown & Bacon, 2009; Kong & Carlyon, 2007; Qin & Oxenham, 2006; Zhang, Dorman, & Spahr, 2010).

### Binaural Hearing Benefits

Preservation of low frequency hearing in the implanted ear also provides potential benefits through enabling provision of acoustic stimulation to both ears. Five out of the 18 inclusion studies examined the impact of utilizing both ears for binaural benefit in speech perception and localization. Three studies that evaluated word performance in quiet (Gantz et al., 2005; Lorens et al., 2008; Simpson et al., 2009) reported higher speech perception scores in the binaural condition (CI+HAi + HAc opposite ear) compared to the monaural

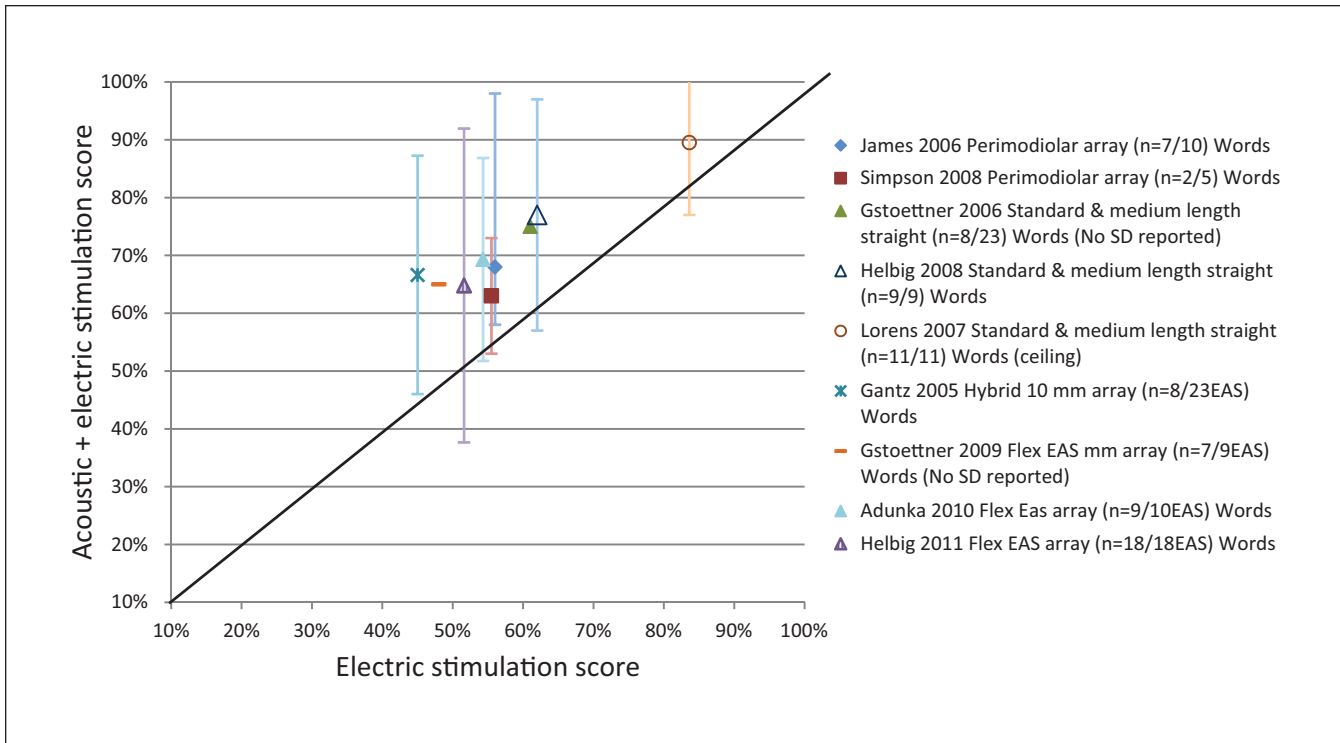

**Figure 3a.** EAS benefit for implanted ear in Quiet: Electric only stimulation (CI) versus Acoustic & Electric Stimulation (CI+HAI) word scores (Mean  $\pm$  SD) for various electrode array studies.

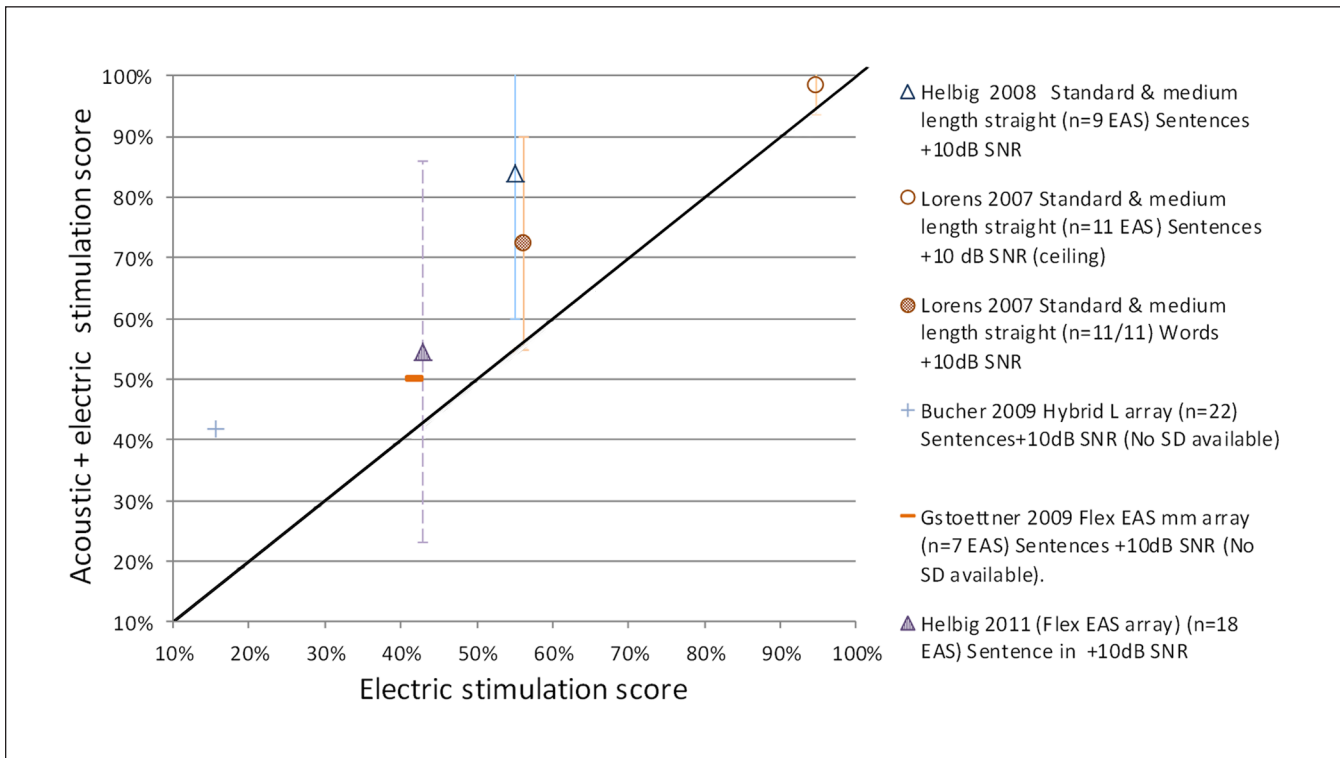

**Figure 3b.** EAS benefit for implanted ear in Noise: Electric only stimulation (CI) versus Acoustic & Electric Stimulation (CI+HAI) scores (Mean  $\pm$  SD) for various electrode array studies.

condition (CI+HAi) on average, however, the studies did not reach statistical significance. Two studies evaluated sentence perception using an adaptive test with spatially separated target speech and noise sources (Dunn et al., 2010; Gifford et al., 2010). Binaural hearing did however provide a statistically significant benefit for sentence perception in spatially separated speech and noise sources. The binaural improvement for speech reception thresholds (SRTs) was about 3dB in both studies. Dunn and colleagues (2010) proposed that the potential advantage lies in the utilization of two similar acoustic signals with identical sound processing being presented to the central auditory system ear. One study to date that reported on horizontal localization found a significant difference in localization performance for the binaural condition compared to the monaural or bimodal condition (Dunn et al., 2010). The author concluded that the bilateral hearing aids providing similar processing from both ears improved the perception of audible low frequency input below 800 Hz, where primarily interaural time level differences (ITD) cues help to localize the sound source. The ability to localize the sounds was dependent on participants being able to perceive sounds in both ears with the acoustic amplification provided by the hearing aids, and independent of whether there was electrical stimulation or not.

Only a small number of studies have utilized test arrangements to evaluate binaural benefit derived through the provision of acoustic stimulation to both ears. The review highlights the need for the inclusion of more studies investigating and using outcomes measures that better quantify binaural advantages such as; spatially separated speech and noise, localization, and also functional benefit in everyday life or Effective (EV) evidence, rather than traditional outcome measures utilized in standard cochlear implantation.

### **Q3. What clinical fitting practices have been developed for devices that combine acoustic and electric stimulation?**

EAS outcome studies report utilizing a range of fitting practices for both the acoustic and electric stimulation. A summary of the 16 studies selected to address the third clinical question are shown in Tables 4.

#### ***Low Frequency Acoustic Amplification for Combined Acoustic and Electric Stimulation***

The provision of acoustic amplification to the residual low frequency hearing preserved in the implanted ear in EAS studies has been either via an independently operating, commercially available digital in-the-ear hearing aid worn together with a behind-the-ear cochlear implant sound processor, or more recently, via EAS sound processors that deliver both acoustic and electrical stimulation through one

integrated device. A wide range of conventional hearing aid prescriptive rules have been used in the studies to determine the amplification applied to the low frequency hearing in the implanted ear. Procedures for linear and nonlinear amplification such as: Half gain rule (Kiefer et al., 2005; Helbig et al., 2011), NAL-RP (Gantz & Turner, 2009; Dunn et al., 2010), DSL [I/O] (James et al., 2005; Fraysse et al., 2006) and NAL-NL1 (Karsten et al., 2012, Simpson et al., 2009) have been applied. The majority of studies apply a variation to these standard hearing aid fitting rationales for the high frequency region of the audiogram. In EAS device fittings the high frequency information is provided via electrical stimulation via the cochlear implant part, therefore, acoustic amplification is only provided up to the point where there is useful residual hearing or cut off where hearing is not measurable. This frequency up to which acoustic amplification is provided ranged from the frequencies where the hearing threshold exceeded 60 dB HL (Helbig et al., 2009; Nopp & Polak, 2010); exceeded 80 dB HL (Fraysse et al., 2006; Gstoettner et al., 2008; James et al., 2005; Lenarz et al., 2009) or 85 dB HL (Gantz et al., 2009; Simpson et al., 2009) on the individual's audiogram. Practical considerations such as prevention of feedback and off-frequency listening, providing sufficient loudness, improving mould or loudness comfort and prolonging battery life, have also been reported as the reason for applying variation to the standard fitting rationale.

Only two studies reported on the broadband amplification applied to the contralateral ear for subjects in their study (Dunn et al., 2010; Gantz et al., 2009). Hearing aid verification using real ear probe measures to obtain objective measures of hearing aid fitting were reported in only two studies (Dunn et al., 2010; Vermeire, Anderson, Flynn, & Van de Heyning, 2008). Four studies report performing loudness balancing tasks between the electric and acoustic signal in the implanted ear. Loudness rating scales were used and adjustments to gain of hearing aid were made in response to individual feedback to live voice recorded speech or noise signals at different input levels (Karsten et al., 2012; Kiefer et al., 2005; Simpson et al., 2009). One study reported using a loudness scaling test performed at 500 Hz with the CI component deactivated to determine the setting of the AGC threshold (Helbig et al., 2011). The majority of the studies examined did not provide any information on whether a loudness balancing procedure between the electric and acoustic signal was used.

To date, two studies have evaluated the influence of different acoustic fitting parameters in the implanted ear on EAS outcomes (Vermeire et al., 2008; Polak, Lorens, Helbig, McDonald, & Vermeire, 2010). A small study by Vermeire and colleagues (2008) investigated the effect of two parameters; acoustic amplification (gain) and frequency range of amplification on speech performance, and subjective



**Table 4.** A Summary of 16 EAS Articles That Provide Information on Fitting Practices for Devices That Combine Electric and Acoustic Stimulation.

Study	Acoustic parameters	Electric parameters	Key findings	Level, grade strength, EV/EF	MD score
1. Perimodiolar electrode array ("Advance off stylet" technique) studies					
Frayse et al., 2006	N = 9 adults. Hearing aid: digital, high-power Phonak ITE Prescription: DSL(i/o) for thresholds ≤ 80 dB HL. Gain reduced > 80dB HL to improve loudness comfort & prevent feedback. Small vent used to remove occlusion discomfort where possible.	Sound Processor: ESPrit 3G, ACE, 20 channels, 8 maxima & Channel stimulation rate: 900 pps. Programs used: 1. Unrestricted electric bandwidth: Full frequency range (188 - 8000 Hz). 2. Restricted electric bandwidth & disabling low frequency channels and shifting frequency-to-electrode allocation. * 1 month each program (randomized) then both available.	Equivalent speech scores for both programs, however, a trend for better scores in noise for restricted electrical bandwidth. Considerable variability in individual performance. Preference: 7/9 preferred restricted electric bandwidth program 2/9 preferred unrestricted electric bandwidth (increasing/ fluctuating HTLs).	3B II EF	27/48
Simpson, McDermott, Dowell, Sucher, & Briggs, 2009	N = 5 adults Hearing aid : digital, high-power Phonak ITE/BTE Prescription: NAL-NL I for thresholds ≤ 90 dB HL. Gain reduced to minimum > 90 dB HL (where hearing considered no longer "useful." Loudness balance task: 9 point loudness rating scale. Subjects listened to running speech at 60dBA to ensure "comfortable" & ICRA noise at 80dB A to ensure "loud but comfortable"	Sound Processor: Freedom, 22channels, ACE, 900Hz, 8 maxima MP stimulation. Programs used: 1. Unrestricted electric bandwidth (250 - 7438Hz) 2. Restricted electric bandwidth Pitch-matched: Pitch-scaling procedure to determine lower frequency boundary for electric stimulation. ABA design (A pitch-matched 12 weeks, B conventional 12 weeks & A pitch matched for 2 weeks)	No significant difference with unrestricted or pitch-mapped programs for group words or SNRs cores (n = 5). Participants using combined electric and acoustic stimulation in the implanted ear (n = 2) obtained best score with restricted electrical bandwidth programs.	3B II EF	30/48
2. Standard length (Reduced/shallow insertion depth) and/or medium length straight electrode array studies					
Kiefer et al., 2005	N = 13 adults Hearing aid: nonlinear compression high-power ITE Prescription: Half gain rule in the 125 to 1000Hz range only. Standard fitting procedures not implemented as it produced "unsatisfactory results where high levels of amplification in mid-high frequency region were either not possible, or desired by participants." Loudness balanced Loudness scaled and individual adjustments made when necessary.	Electric parameters Sound Processor: TEMPO+, CIS strategy, 7-11 active channels available. Total stimulation rate: 18180 pps Programs used: 1. Unrestricted electric bandwidth (300 - 5500 Hz). 2. Restricted electric bandwidth (650 Hz - 5500Hz) 3. Restricted electric bandwidth (1000 Hz - 5500Hz) Initially fit with unrestricted BW program in CI alone condition until stable thresholds reached. Then all 3 programs used for 2-3 weeks prior to testing.	12/13 participants had best results with & chose the full frequency unrestricted electric bandwidth program (300 - 5500Hz) 1/13 obtained best score with restricted program (650Hz -5.500Hz) Authors comment : "Overlap not detrimental" Participants received longer exposure to full frequency range program versus others	3C II EF	24/48

(continued)

**Table 4. (continued)**

2. Standard length (Reduced/shallow insertion depth) and/or medium length straight electrode array studies		Level, grade strength, EV/EF	MD score
Study	Acoustic parameters	Electric parameters	Key findings
Skarzynski, Lorens, Piotrowska, & Anderson, 2007	N = 9 children used residual hearing with CI. 4/9 used a DUET combined EAS stimulation device. 1/9 uses CI only. Pediatric Fitting: "The hearing aid component was fitted according to the required gain, using the half-gain rule (where the threshold of a given frequency is divided in half to obtain the predicted value). The recommended fitting procedures were followed and no special consideration had to be made for fitting the children."	Sound Processor: Duet EAS (which uses a TEMPO+ processor and a 2 channel hearing aid in one device.) 12 channels' Programs used: 1. Restricted electric bandwidth between 500 & 1000Hz. Frequency cutoff point determined by audiogram. (n = 8) 2. Unrestricted electric bandwidth (n = 1) range not reported	Author's comments: "The ci was programmed in such a way that there is no overlap with acoustic perception, so as to not interfere with this perception." "The children readily accepted the combination of electric and acoustic information. These results suggest that there are no special issues when fitting children their natural acoustic low frequency hearing"
Gstoettner et al., 2008	N = 18 adults Hearing aid: digital high-powered Oticon ITE: 2 channel, seven band instrument & peak gain 60dB & MPO 127 dB SPL Prescription: Not reported. Amplification provided for all low frequencies where audiogram ≤ 80 dB HL. For > 80dB HL input into the software arbitrary 10dB HL threshold * Aids fitted 2 months post-activation.	Sound Processor: TEMPO+, CIS strategy, Total stimulation rate: 18180 pps Programs used: 1. Unrestricted electric bandwidth (up to 7000 Hz) for 2 months. 2. Restricted electric bandwidth for Cross-over frequency determined by the frequency at which threshold >65 dB HL on unaided audiogram. Initially used unrestricted BW program in CI alone condition for 2 months & then restricted program fit.	12/18 participants had sufficient residual hearing to utilize ipsilateral electric and acoustic stimulation. Author's comments: "Six however felt that additional acoustic input was either detrimental to their sound quality or insufficient in terms of amplification, or they disliked management of the two devices and thus choosing not to use EAS."
Vermeire, Anderson, Flynn, & Van de Heyning, 2008	N = 4 adults Hearing aid: digital high-power Oticon ITE: Vent as large as possible to reduce occlusion effect. Prescription: Half gain rule & VoiceFinder algorithm enabled which adjusts gain based on input signal. HA Programs used: 1. For frequencies ≥ 85dB HL a 10 dB HL threshold was input into audiogram (frequencies ≥ 85dB HL not amplified) 2. For frequencies ≥ 120dB H a 10 dB HL threshold was input into audiogram (frequencies ≤ 115 Hz amplified). 3. Gain & compression predicted by Adapto fast fitting software. 4. Gain increased by 6db SPL & compression reduced for a LF 'boost.'	Sound Processor: TEMPO+ sound processor Programs used: 1. Unrestricted electric bandwidth stimulation (200 - 7000Hz) 2. Restricted electric bandwidth Cross-over frequency determined by the frequency at which threshold ≥ 65 dB HL on unaided audiogram Fitting: Participants initially fit with unrestricted BW program for 2 months. Then 30 minute adjustment period with each listening program & 1 practice list for each section.	Acoustic signal: Highest level of performance in noise & subjective reports of listening effort (Visual Analogue Scale) when frequencies ≤ 120 dB HL amplified & 6dB gain (low boost). Electric signal: "Best results" in noise (+15 & +10 SNR) with restricted electrical bandwidth program. 3/4 subjects with steeply sloping losses performed better with restricted electric bandwidth program & low boost to gain. One subject with flat severe loss performed better with unrestricted program. NSD effects for any condition on VAS scale.

(continued)

**Table 4. (continued)**

2. Standard length (Reduced/shallow insertion depth) and/or medium length straight electrode array studies		Level, grade strength, EV/EF	MD score
Study	Acoustic parameters	Electric parameters	Key findings
Helbig, Baumann, Helbig, von Malsen-Waldkirch, & Gstoettner, 2008	<p>N = 9 adults</p> <p>Acoustic via ITE and/or DUET EAS.</p> <p>Prescription: Half gain rule &amp; optimized according to subject's responses.</p> <p>Fitting Issues: "The handling of two separate devices with different batteries and with different battery life spans is cumbersome and was reason enough to reject the use of the hearing aid. In patients with poorer low frequency thresholds, the ITE amplification was often insufficient, especially in the frequencies below 500Hz, which made fitting to provide sufficient gain very difficult."</p>	<p>Sound Processor: Duet EAS</p> <p>Programs used: Not reported</p>	<p>Initially, 9 participants wore ITE but then 6 rejected it.</p> <p>Author's comment: "Despite the fact that an increase in speech acoustic amplification, this lack of gain and wearing difficulties resulted in individuals discontinuing to use the hearing aid before they could possibly perceive its benefit."</p>
Helbig & Baumann, 2009	<p>N = 15 adults</p> <p>Acoustic: DUET 2 channel hearing aid with 4 trimming controllers: gain, LFS, volume and AGC. Gain at 500 Hz is adjustable between 27-42 dB HL. Peak gain is 52 dB HL at 1.6KHz.</p> <p>Prescription: Half gain rule set based on hearing threshold at 500Hz,</p> <p>Slope: "Low frequency slope (LFS) set between 0-18dB/octave. Subtraction of thresholds at 250 &amp; 500 Hz. Half of this difference determines LFS.</p> <p>Loudness balance: volume controller used for overall attenuation &amp; ranges between 0-47 dB HL.</p> <p>Loudness scaling test performed at 500Hz (acoustic only) to determine setting of AGC threshold (input levels between 40-70 dB HL)</p>	<p>Sound Processor: Duet EAS</p> <p>Program used:</p> <p>1. Restricted electric bandwidth</p> <p>Cross-over frequency determined by the frequency at which threshold is &gt;65 dB HL on unaided audiogram.</p> <p>Fitting issues: Larger overlap of acoustic &amp; electric "individuals with nearly normal hg in low frequency complained of "echo."</p> <p>"Amplification peak at 1.6KHz led to problems of "distorted speech" dead regions in the area and lower regions activated."</p>	<p>11 of 15 subjects continued to use the combined DUET EAS sp .</p> <p>4 of 15 participants rejected acoustic amplification due to limited benefit.</p> <p>"The "DUET user" group had residual hearing that was better than 75dB in the 500Hz frequency region or below.</p> <p>The "non-user" group had a hearing loss of greater than 55dB HL at 125 Hz, 70 dB HL at 250 Hz and 98 dB HL at 500Hz</p> <p>Acceptance of acoustic amplification was dependent on postoperative hearing thresholds"</p>
Polak, Lorenz, Helbig, McDonald, S., & Vermeire, 2010	<p>N = 24 adults</p> <p>Acoustic 2 channel hearing aid via DUET™ EAS processor</p> <p>Prescription: Initially Half gain rule plus individual adaptation.</p> <p>3 Acoustic Parameters varied:</p> <p>LF slope: <math>Th_{500} / Th_{250} / o, 18dB/octave</math></p> <p>Compression threshold: 40.55,70 dB</p> <p>Compression: ratio: 1:1, 1:1.33, 1:2</p>	<p>Sound Processor: Duet EAS</p> <p>Programs used:</p> <p>All restricted electric bandwidth programs: 200Hz from unaided audiogram at:</p> <ol style="list-style-type: none"> <li>≥ 50 dB HL threshold</li> <li>≥ 65 dB HL threshold</li> <li>≥ 80 dB HL threshold</li> </ol> <p>Fitting: Participants had 2 hr or 1 day adjustment period with each listening program.</p>	<p>Study reported optimization of the parameters had a 'strong effect' on speech performance, however, no statistical analysis of the results.</p> <p>Parameters that had most impact on "overall" benefit were compression ratio &amp; frequency where electrical stimulation started. No details on which of the parameters varied were optimal</p>

(continued)

**Table 4. (continued)**

Study	Acoustic parameters	Electric parameters	Key findings	Level, grade strength, EV/EF	MD score
2. Standard length (Reduced/shallow insertion depth) and/or medium length straight electrode array studies					
Nopp & Polak, 2010	Review paper reports on Polak study (2007) N = 15 adults Acoustic via DUET™ EAS processor Prescription: Half gain rule in the region of 125 to 1500Hz Hz, plus individual adaptation (slope, AGC threshold, volume).	Sound Processor: Duet EAS Programs used: 1. Unrestricted electric bandwidth: (range 200 - 8.5000 Hz). Restricted electric bandwidths at: 2. ≥ 50 dB HL threshold 3. ≥ 65 dB HL threshold 4. ≥ 80 dB HL threshold Fitting: 1 day period for adjustment with each program	On average, best speech understanding achieved when restricted electric bandwidth range set to the frequency at 65HL on the unaided audiogram. 78% subjects had poorer speech test scores with unrestricted electric bandwidth (full frequency range) program. No data provided	6D III Review paper	*
Skarzynski & Lorens, 2010	Design: N = 25 total 2 children groups: 1. Partial Deafness (n = 15) 2. Platinum HA users (n = 10) Acoustic: Different devices used*	Sound Processor: Duet EAS with 8-11 active channels Program used: 1. Restricted electric bandwidth. Cross-over frequency determined by unaided audiogram. Low frequency ranged from 300 & 1000 Hz. The upper frequency was 8.5 kHz. Fitting: No extra-cochlea electrodes used. Determined by impedance telemetry & hearing sensation reports.	Pediatric Fitting: "slight overlap with acoustic perception" * Individuals used either their natural acoustic low frequency hearing together with CI for everyday listening or used a DUET to amplify low frequencies and ci to amplify high frequencies in same ear (n = 22/25). Participants with "nonfunctional partial preservation" used ci only (n = 3/25).	4C II EF	25/48
3. Straight electroacoustic electrode array studies					
Skarzynski, Lorens, Piotrowska, & Skarzynski, 2010	Design: N = 95 total 2 group: Adults (n = 62) Children (n = 32) Acoustic: Different devices used Individuals used either their natural acoustic low frequency hearing together with CI or used a DUET to amplify low frequencies and ci to amplify high frequencies in same ear or CI alone.	Sound Processor: Duet EAS Programs: Not reported "The Skarzynski PDT Classification. A proposed classification system to assist with decision making process for amplification in the IE based on individual's postoperative PTA at 125 Hz, 250 Hz & 500Hz.	3 Groups: A. Electrical component (EC) & nonamplified residual hearing for patients with normal or slightly elevated hearing losses (n = 30). B. Electric-acoustic stimulation (EAS) for patients with mild-to-severe loss (n = 43). C. Electrical component (EC) only for patients with nonfunctional preservation or complete hg loss(n = 22).	4C II EF	20/48

(continued)

Table 4. (continued)

3. Straight electroacoustic electrode array studies		Level, grade strength, EV/EF	MD score
Study	Acoustic parameters	Electric parameters	Key findings
Gantz et al., 2009	N = 87 adults Acoustic: Some used residual hearing & some used different devices: Prescription: NAL, Amplification provided up to a cutoff of useful residual hearing $\leq 90$ dB HL. Broadband amplification applied to the contralateral ear for all subjects (NAL). 2 subjects preferred to use their residual low frequency hearing in both ears rather than hearing aids (both had mild to moderate hearing thresholds up to 750Hz in both ear). 2 used amplification only in contralateral ear (1 profound loss & 1 with complete loss.)	Sound Processor: All processors programmed with a standard CIS type of processing strategy, 6 channels available. Program used: 1. Restricted electric bandwidth Cross-over frequency determined by frequency at which the threshold is $\geq 90$ dB HL on unaided audiogram. Frequency range of 688 - 7938 or 1063 - 7938 Hz	3C II EF 24/48 Author comments: "In cases where a subject lost sufficient hearing so that amplification was not beneficial, the implant was programmed to correspond with the upper frequency cutoff of the contralateral ear as well as a broadband MAP." Testing proceeded with the patient's preferred map, which in all cases was the map with the abbreviated frequency allocation.
Lenarz et al., 2009	N = 32 adults Acoustic: ITE in IE: Amplification provided for all low frequencies up to the frequency point where audiogram was 80 dB HL. Recipients not using the HA were fitted with a standard CI program covering the low frequency range.	Sound Processor: Freedom. 22 channels available Program used: 1. Restricted electric bandwidth. Cross-over frequency determined by frequency at which the threshold is $\geq 80$ dB HL on unaided audiogram.	3C II EF 28/48 25 participants used acoustic and electric stimulation in the implanted ear. 6 used CI alone in the implanted ear. 1 preferred to used residual hearing without a HA in the implanted ear. 8/32 participants reported basal most electrodes cause non-auditory sensations (pain or facial nerve stimulation) All sensations were resolved when electrodes deactivated.
Dunn, Perreau, Gantz, & Tyler, 2010	N = 11 adults Hearing aids: I5 channel Phonak ITE in ipsilateral ear and own BTE in contralateral ear to implant. Prescription: NAL-RP targets at all low frequencies bilaterally. Fitting: Real ear probe measurements to verify targets matched.	Sound Processor: Stimulation rates of up to 2400 Hz or 3500 Hz pulses per channel. 6 channels available. Program used: 1. Restricted electric bandwidth "Cochlear implant frequency response was set to supplement the subject's acoustic hearing as determined by their audiogram."	3C II EF 28/48 Author comments: "It appears when listeners have similar processing bilaterally through the use of bilateral hearing aids, listeners are able to take advantage of ITD cues, and the addition of the cochlear implant did not disrupt overall performance."
Karsten et al., 2012	N = 10 adults. Prescription: NAL-NLI targets Fitting: Real-ear aided response measures to verify targets matched. Loudness balancing using live speech	Sound Processor: Nucleus Hybrid Stimulation rates of 900 to 2500 Hz per channel, ACE processing strategy, 6-10 channels available. Program used: Restricted electric bandwidths: Cross-over frequency determined by frequency at which Electric and acoustic signal: 1. Meet 2. Overlap 3. Gap between Fitting: Minimum of 11 days for adjustment with each program	3B II EF 30/48 Significantly better sentence scores in noise & subjective ratings for the "meet" non-overlapping and restricted bandwidth program. No significant difference between programs for consonant recognition in quiet.

## Abbreviations:

ACE= Advanced combination encoder; BTE= Behind the ear; CI: Cochlear implant; dB = decibel; CIS= Continuous interleaved sampling DSL = Desired Sensation Level; DSLi/o = Desired Sensation Level Input/Output; EF=Effacious or laboratory based measures; EV= Effectiveness or real world measures; HA = hearing aids; IE= Implanted ear; ITE=In the ear; HA= Hearing aids; MD= MacDermid score;MP= Monopolar stimulation; NAL RP = National Acoustic Laboratories-Revised Profound; NAL-NLI = National Acoustic Laboratories-Non-linear I; NIE = nonimplanted ear; NH= Normal hearing; PDCI=Partial deafness cochlear implantation; PTA= Pure tone Audiometry; PTT=Pure tone thresholds; SNR= Sensorineural hearing loss; SNR = signal to noise ratio; PPS=Pulses per second.

measure of benefit in 4 adults with a MED-EL M- electrode array and ITE hearing aid in the same ear. Results showed that additional low frequency gain (+6dB SPL in the low frequency range) and providing a wider frequency range for amplification (amplification up to thresholds of 120dB HL on the audiogram) provided the highest mean speech scores. Unfortunately, small participant numbers and limited listening experience of 30 min with each condition limit the strength of evidence and therefore the ability to support definitive conclusions. A second study by Polak et al. (2010) evaluated compression, compression threshold, and low frequency slope in the implanted ear (nonimplanted ear plugged) in 24 EAS participants. The authors reported that whilst all the parameters tested had an effect on speech perception in quiet and noise in the implanted ear, the parameters that had the most impact on overall benefit was the compression ratio and the frequency where electrical stimulation started. No statistical analysis of the results was reported to evaluate significance. Limited listening experience with each condition (2 hr or 1 day) also limits the strength of evidence from this study. Further investigation is required. A summary of hearing aid parameters used in EAS fittings is shown in Table 3.

### *Discontinuation of Low Frequency Acoustic Amplification Postoperatively*

A number of studies report discontinuation or rejection of amplification in the implanted ear by some individuals. Typically, they can be divided into three groups. The first group consisted of individuals with a severe to profound hearing loss in the implanted ear (typically >80 dB HL across the frequency range) (Gantz et al., 2009; Gstoettner et al., 2008; Helbig et al., 2011; Simpson et al., 2009; Skarzynski et al., 2010) and without sufficient useful hearing for combined electric acoustic stimulation. It is not clear whether this is due to limited benefit or the “perceived” limited benefit by the recipient or as a result of practical issues concerning physical or loudness discomfort around fitting such losses, the inconvenience of wearing an additional component or a combination of factors. In addition Helbig et al. (2011) reported that “amplification was often reported to be insufficient, especially in the frequencies below 500Hz, which made fitting to provide sufficient gain very difficult.” A study by Helbig and Baumann (2009) examined the relationship between the postoperative pure tone audiometric thresholds in the implanted ear and wearer acceptance of a combined DUET EAS sound processor in 15 individuals. The study found that the 11 participants continued to use combined electric and acoustic stimulation, had residual hearing that was <75dB HL in the 500Hz frequency region or below. The remaining four participants who rejected acoustic amplification, had a hearing loss of greater than 55dB HL at 125 Hz, 70 dB HL at 250 Hz and 98 dB HL at 500Hz. The authors deduced that acceptance of acoustic amplification was dependent on postoperative hearing thresholds.

The second group that did not use amplification in the implanted ear consists of individual with significant amounts of residual hearing between 125-500Hz and PTA within the normal to mild range (Gantz et al., 2009; Lenarz et al., 2009; Skarzynski, Lorens, Piotrowska, & Anderson, 2007; Skarzynski & Lorens, 2010). This group did not require amplification as they could combine their natural acoustic hearing with electric stimulation. Skarzynski et al. (2010) proposed a classification system where electrical stimulation alone should be used in the implanted ear if low frequency hearing (125 Hz - 500 Hz) is  $\leq 30$  dB HL or  $\geq 80$  dB HL, and acoustic amplification considered in the implanted ear for all others.

The final group consists of individuals that report management and wearing of a supplementary acoustic instrument was cumbersome, inconvenient or contraindicated due to issues of the external ear canal such as acute external otitis media and was the reason to reject the use of the hearing aid (Gstoettner et al., 2008; Helbig, Baumann, Helbig, von Malsen-Waldkirch, & Gstoettner, 2008; Helbig et al., 2011) despite some having sufficient residual hearing to utilize and potentially benefit from ipsilateral electric and acoustic stimulation. Whilst the introduction of less cumbersome EAS sound processors which deliver both acoustic and electrical stimulation in one integrated device has helped reduce some of these wearing issues, there are still a number of participants that have discontinued amplification (Helbig et al., 2008; Helbig & Baumann, 2009). Further systematic investigation of fitting procedures aimed at maximizing outcomes for recipients of electric-acoustic devices could potentially facilitate the continued use of acoustic amplification following implantation, reduce individual variability with regard EAS benefit, increase recipient outcomes and wearer satisfaction.

### *Cochlear Implant Electric Fitting for Combined Acoustic and Electric Stimulation*

To date, the electric parameters selected for electric stimulation in EAS applications are generally default settings determined by the cochlear implant system used. The default parameters utilized in studies include sound coding strategies, stimulation rates, number of channels. A summary of electric parameters used in EAS fittings is shown in Table 3.

One electric fitting parameter that has received a considerable amount of attention is the electric-to-acoustic frequency boundary allocation, also known as the crossover frequency. The crossover frequency refers to the way in which frequencies in the input signal are selected and mapped for acoustic stimulation and for electric stimulation in the implanted ear (the frequency at which stimulation changes from acoustic to electrical). There is potential to provide overlapping inputs for low frequencies with the CI sound processor via electrical stimulation (200 Hz - 1000Hz) and the residual low frequency acoustic hearing

(0 Hz -1000 Hz). Therefore, the same acoustic frequency can be represented in different cochlea locations via two very different modalities within the implanted ear. A number of studies attempted to quantify how much overlapping, if any, of the frequency range would provide the best outcomes for speech perception. The majority of published evidence that suggests that limiting or completely restricting the range of incoming stimulus frequencies that could be represented both acoustically and electrically provided potential speech perception benefits and was subjectively preferred by the majority of listeners. There was also a trend for better speech perception scores in noise when a restricted electrical bandwidth was used in the CI program. (Frayse et al., 2006; Helbig et al., 2011; James et al., 2005; Karsten et al., 2012; Nopp & Polak, 2010; Simpson et al., 2009; Vermeire et al., 2008). The minority of listeners reported in the studies that preferred full electrical frequency ranges and overlapping acoustic and electric inputs, typically had more severe to profound losses or fluctuating losses compared to the steeply sloping losses (Frayse, et al., 2006; Gantz et al., 2009; Vermeire et al., 2008; Nopp & Polak, 2010). While listeners with normal to mild hearing thresholds in the low frequency region reported “echo” and “distorted” sound quality with full unrestricted electric or overlapping acoustic and electric programs in a few studies (Helbig & Baumann, 2009; Karsten et al., 2012).

The general approach for fitting these devices has been based on the principle that electric stimulation by the CI processor should begin at the point where acoustic hearing is no longer beneficial. The method used to determine this input frequency range however, has varied widely in studies reviewed. The majority of studies report crossover frequency selection being based on the individual’s unaided audiogram and where the dead regions of the cochlea might limit benefit. In these studies the selection of the crossover frequency based on the audiogram can be grouped into two: the frequency at which the individual’s thresholds are  $\geq 65$  dB HL on the audiogram (Helbig et al., 2011; Nopp & Polak, 2010) or  $\geq 80$  dB HL on the audiogram (Frayse et al., 2006; Gantz et al., 2009; Gstoettner et al., 2008; James et al., 2005; Lenarz et al., 2009; Simpson et al., 2009). In addition, Greenwood’s frequency-to-place map calculations (Kiefer et al., 2005) and individually optimized frequency derived from psychophysical measurements (Baumann, Rader, Helbig, & Bahmer, 2011; Simpson et al., 2009) methods have also been reported. The findings do not provide a clear understanding of where to set the restricted electric bandwidth range based on an individual’s hearing. All studies have used sentence or monosyllabic word recognition in quiet and noise as outcome measurements. Although several studies include comments on participant’s subjective preference, only two study (Karsten et al., 2012; Vermeire et al., 2008) used a qualitative rating scales such as the visual analogue scale (VAS) to rate “ease of listening” effort during testing or in everyday life. The

question of whether systematic adjustment of electric and acoustic output levels and parameters improved music perception, localization or functional performance in real life usage outcome, remains unclear.

### *Pediatric Fitting for Combined Electric and Acoustic Stimulation*

Finally, two studies by Skarzynski and colleagues (2007), and Skarzynski & Lorens (2010) describe fitting procedures for EAS devices in the pediatric population. Skarzynski et al. (2007) described using a standard fitting procedure developed for adults when fitting children. No special consideration had to be made and there were no specific issues relating to these fittings. The author stated that “children readily accepted the combination of electric and acoustic information.”

Whilst there is some evidence that suggests a systematic fine-tuning of the electric and acoustic stimulation may be beneficial, further research is needed to investigate fitting procedures to maximize outcomes for adults and children with different degrees of residual hearing who use EAS implant systems.

### **Conclusion**

The preservation of low frequency residual hearing following cochlear implantation provides the capacity for combining electric and acoustic stimulation in the implanted ear to benefit the user. This emphasizes the importance of continued development of electrode design and surgical techniques aimed at preserving residual low frequency hearing in the implanted ear. This knowledge also highlights the need for systematic reporting procedure for all postoperative hearing preservation studies and underscores the importance of appropriate preoperative counseling to ensure realistic expectations of potential benefits versus and risks.

The EAS research evidence to date has also shown that benefit from combining electric and acoustic stimulation in both the monaural and binaural condition can be consistently demonstrated over electric alone condition. This potential EAS advantage highlights the value of the endeavors to preserve low frequency hearing in the implanted ear and to continue amplification where appropriate. Consideration for the inclusion of additional outcome measures that help effectively evaluate the benefit of the additional acoustic low frequency hearing in both ears for speech perception, localization and functional performance everyday life is crucial.

A wide variety of acoustic and electric fitting parameters and methods have been used and the individual variations applied to amplification of the high frequencies. Further research is also needed to investigate fitting procedures to maximize outcomes for recipients with different degrees of hearing loss who use devices that combine electric and

acoustic stimulation. Ensuring that the electric and acoustic stimulation devices are fitted optimally through empirically validated procedures and incorporating current knowledge of potential binaural benefit into a comprehensive individual rehabilitation program will facilitate their continued use following implantation and maximize outcomes.

### Acknowledgment

The authors are grateful for the financial support of the Commonwealth of Australia through the HEARing Cooperative Research Centre, an Australian Government initiative.

### Declaration of Conflicting Interests

The authors declared no potential conflicts of interest with respect to the authorship and publication of this article.

### Funding

The authors disclosed receipt of the following financial support for the research and authorship of this article: This work is supported by the Commonwealth of Australia through the HEARing Cooperative Research Centre.

### References

- Adunka, O. F., Pillsbury, H. C., Adunka, M. C., & Buchman, C. A. (2010). Is electric acoustic stimulation better than conventional cochlear implantation for speech perception in quiet? *Otology & Neurotology*, *31*, 1049-1054
- Arnoldner, C., Helbig, S., Wagenblast, J., Baumgartner, W., Hamzavi, J., Riss, D., & Gstöettner, W. (2010). Electric acoustic stimulation in patients with postlingual severe high-frequency hearing loss: Clinical experience. *Advances in Oto-Rhino-Laryngology*, *67*, 116-124.
- Baumann, U., Rader, T., Helbig, S., & Bahmer, A. (2011). Pitch matching psychometrics in electric acoustic stimulation. *Ear & Hearing*, *32*, 656-662.
- Brockmeier, S. J., Peterreins, M., Lorens, A., Vermeire, K., Helbig, S., Anderson, I., & Kiefer, J. (2010). Music perception in electric acoustic stimulation users as assessed by the Mu.S.I.C. test. *Advances In Oto-Rhino-Laryngology*, *67*, 70-80.
- Brown, C. A., & Bacon, S. P. (2009). Achieving electric-acoustic benefit with a modulated tone. *Ear Hear*, *30*, 489-493.
- Büchner, A., Schüssler, M., Battmer, R. D., Stöver, T., Lesinski-Schiedat, A., & Lenarz, T. (2009). Impact of low-frequency hearing. *Audiology and Neurotology*, *14*(Supplement 1), 8-13.
- Cox, R. M. (2005). Evidence-based practice in provision of amplification. *Journal of the American Academy of Audiology*, *16*, 419-438.
- Dorman, M. F., Gifford, R. H., Spahr, A. J., & McKarns, S. A. (2008). The benefits of combining acoustic and electric stimulation for the recognition of speech, voice and melodies. *Audiology and Neurotology*, *13*(2), 105-112.
- Dorman, M. F., & Gifford, R. H. (2010). Combining acoustic and electric stimulation in the service of speech recognition. *International Journal of Audiology*, *49*, 912-919.
- Driver, S., & Stark, E. (2010). Real life outcomes for adult EAS patients and adult traditional CI patients. *Cochlear Implants International*, *11*(3), 113-118.
- Dunn, C. C., Perreau, A., Gantz, B., & Tyler, R. S. (2010). Benefits of localization and speech perception with multiple noise sources in listeners with a short-electrode cochlear implant. *Journal of the American Academy of Audiology*, *21*, 44-51.
- Fraysse, B., MacÅas, Å. R., Sterkers, O., Burdo, S., Ramsden, R., Deguine, O., & James, C. (2006). Residual hearing conservation and electroacoustic stimulation with the Nucleus 24 Contour Advance cochlear implant. *Otology & Neurotology*, *27*, 624-633.
- Friedland, D. R., & Runge-Samuelson, C. (2009). Soft cochlear implantation: Rationale for the surgical approach. *Trends in Amplification*, *13*(2), 124-138.
- Gantz, B. J., & Turner, C. W. (2003). Combining acoustic and electrical hearing. *Laryngoscope*, *113*, 1726-1730.
- Gantz, B., Turner, C., Gfeller, K., & Lowder, M. (2005). Preservation of hearing in cochlear implant surgery: Advantages of combined electrical and acoustical speech processing. *Laryngoscope*, *115*, 796-802.
- Gantz, B. J., Hansen, M. R., Turner, C. W., Oleson, J. J., Reiss, L. A., & Parkinson, A. J. (2009). Hybrid 10 Clinical Trial. *Audiology and Neurotology*, *14*(Supplement 1), 32-38.
- Gfeller, K. C., Olszewski, C., Turner, C., Gantz, B., & Oleson, J. (2006). Music perception with cochlear implants and residual hearing. *Audiology & Neurotology*, *11*(1), 12-15.
- Gfeller, K., Turner, C., Oleson, J., Zhang, Z., Gantz, B., Froman, R., & Olszewski, C. (2007). Accuracy of cochlear implant recipients on pitch perception, melody recognition, and speech reception in noise. *Ear and Hearing*, *28*, 412-423.
- Gifford, R. H., Dorman, M. F., & Brown, C. A. (2010). Psychophysical properties of low-frequency hearing: Implications for perceiving speech and music via electric and acoustic stimulation. *Advances in Oto-rhino-laryngology*, *67*, 51-60.
- Gstöettner, W. K., Helbig, S., Maier, N., Kiefer, J., Radeloff, A., & Adunka, O. F. (2006). Ipsilateral electric acoustic stimulation of the auditory system: Results of long-term hearing preservation. *Audiology & Neuro-otology*, *1*, 49-56.
- Gstöettner, W. K., Van de Heyning, P., Fitzgerald O'Connor, A., Morera, C., Sainz, M., Vermeire, K., & Adunka, O. (2008). Electric acoustic stimulation of the auditory system: Results of a multi-centre investigation. *Acta Oto-laryngologica*, *128*, 968-975.
- Gstöettner, W. K., Helbig, S., Settevendemie, C., Baumann, U., Wagenblast, J., & Arnoldner, C. (2009). A new electrode for residual hearing preservation in cochlear implantation: First clinical results. *Acta Oto-laryngologica*, *129*, 372-379.
- Gstöettner, W. K., Van de Heyning, P., Fitzgerald O'Connor, A., Kiefer, J., Morera, C., Sainz, M., & Helbig, S. (2011). Assessment of the subjective benefit of electric acoustic stimulation with the abbreviated profile of hearing aid benefit. *Journal for Oto-rhino-laryngology and Its Related Specialities*, *73*, 321-329.
- Helbig, S., & Baumann, U. (2009). Acceptance and fitting of the DUET Device—A combined speech processor for electric acoustic stimulation. *Advances in Oto-rhino-laryngology*, *67*, 81-87.
- Helbig, S., Baumann, U., Helbig, M., von Malsen-Waldkirch, N., & Gstöettner, W. (2008). A new combined speech processor for electric and acoustic stimulation—Eight months experience. *Journal for Oto-rhino-laryngology and Its Related Specialities*, *70*, 359-365.



- Helbig, S., Van de Heyning, P., Kiefer, J., Baumann, U., Kleine-Punte, A., Brockmeier, H., & Gstoettner, W. (2011). Combined electric acoustic stimulation with the PULSARCI(100) implant system using the FLEX(EAS) electrode array. *Acta Oto-laryngologica*, 1-11.
- James, C. J., Albegger, K., Battmer, R., Burdo, S., Deggouj, N., Deguine, O., & Fraysse, B. (2005). Preservation of residual hearing with cochlear implantation: How and why. *Acta Oto-laryngologica*, 125, 481-491.
- James, C. J., Fraysse, B., Deguine, O., Lenarz, T., Mawman, D., Ramos, A., & Sterkers, O. (2006). Combined electroacoustic stimulation in conventional candidates for cochlear implantation. *Audiology and Neurotology*, 11(Suppl. 1), 57-62.
- Karsten, S. A., Turner, C. W., Brown, C. J., Joen, E. K., Abbas, J., & Gantz, B. J. (2012). Optimising the combination of acoustic and electric hearing in the implanted ear. *Ear and Hearing*, Advance Online Publication, doi:10.1097/AUD.0b013e31826a1c99
- Kiefer, J., Pok, M., Adunka, O., Stürzebecher, E., Baumgartner, W., Schmidt, M., & Gstoettner, W. (2005). Combined electric and acoustic stimulation of the auditory system: Results of a clinical study. *Audiology and Neurotology*, 10(3), 134-144.
- Kong, Y. Y., & Carlyon, R. P. (2007). Improved speech recognition in noise in simulated binaurally combined acoustic and electric stimulation. *Journal of the Acoustical Society of America*, 121, 3717-3727
- Lenarz, T., Stöver, T., Buechner, A., Paasche, G., Briggs, R., Risi, F., & Battmer, R. (2006). Temporal bone results and hearing preservation with a new straight electrode. *Audiology and Neurotology*, 11(Supplement 1), 34-41.
- Lenarz, T., Stöver, T., Buechner, A., Lesinski-Schiedat, A., Patrick, J., & Pesch, J. (2009). Hearing conservation surgery using the Hybrid-L Electrode. *Audiology and Neurotology*, 14(Supplement 1), 22-31.
- Lorens, A., Polak, M., Piotrowska, P., & Skarzynski, H. (2008). Outcomes of treatment of partial deafness with cochlear implantation: A DUET Study. *The Laryngoscope*, 118, 288-294.
- MacDermid, J. C. (2004). An introduction to evidence-based practice for hand therapists. *Journal of Hand Therapy*, 17(2), 105-117.
- Mukherjee, P., Uzun-Coruhlu, H., Wong, C. C., Curthoys, I. S., Jones, A. S., & Gibson, W. P. (2012). Assessment of intracochlear trauma caused by the insertion of a new straight research array. *Cochlear Implants International*, 13(3), 156-162.
- Nopp, P., & Polak, M. (2010). From electric acoustic stimulation to improved sound coding in cochlear implants. *Advances in Oto-Rhino-Laryngology*, 67, 88-95.
- Olson, A. D., & Shinn, J. B. (2008). A systematic review to determine the effectiveness of using amplification in conjunction with cochlear implantation. *Journal of the American Academy of Audiology*, 19, 657-671.
- Podskarbi-Fayette, R., Pilka, A., & Skarzynski, H. (2010). Electric stimulation complements functional residual hearing in partial deafness. *Acta Oto-laryngologica*, 130, 888-896.
- Polak, M., Lorens, A., Helbig, S., McDonald, S., & Vermeire, K. (2010). Fitting of the hearing system affects partial deafness cochlear implant performance. *Cochlear Implants International*, 11(Supplement 1), 117-121.
- Qin, M. K., & Oxenham, A. J. (2006). Effects of introducing unprocessed low-frequency information on the reception of envelope-vocoder processed speech. *Journal of the Acoustical Society of America*, 119, 2417-2426.
- Simpson, A., McDermott, H. J., Dowell, R. C., Sucher, C., & Briggs, R. J. S. (2009). Comparison of two frequency-to-electrode maps for acoustic-electric stimulation. *International Journal of Audiology*, 48(2), 63-73.
- Skarzynski, H., Lorens, A., Piotrowska, A., & Anderson, I. (2007). Preservation of low frequency hearing in partial deafness cochlear implantation (PDCI) using the round window surgical approach. *Acta Oto-laryngologica*, 127(1), 41-48.
- Skarzynski, H., Lorens, A., Piotrowska, A., & Podskarbi-Fayette, R. (2009). Results of partial deafness cochlear implantation using various electrode designs. *Audiology and Neurotology*, 14(Supplement 1), 39-45.
- Skarzynski, H. L., & Lorens, A. (2010). Electric acoustic stimulation in children. *Advances in Oto-Rhino-Laryngology*, 67, 135-143.
- Skarzynski, H., Lorens, A., Piotrowska, A., & Skarzynski, P. (2010). Hearing preservation in partial deafness treatment. *Medical Science Monitor*, 16(11), 55-62.
- Skarzynski, H., Lorens, A., Matusiak, M., Porowski, M., Skarzynski, P., & James, C. (2012). Partial deafness treatment with the Nucleus Straight Research Array Cochlear Implant. *Audiology & Neurotology*, 17, 82-91.
- Talbot, K. N., & Hartley, D. E. (2008). Combined electro-acoustic stimulation: A beneficial union? *Clinical Otolaryngology*, 33, 536-545.
- Turner, C. W., Gantz, B. J., Karsten, S., Fowler, J., & Reiss, L. A. (2010). Impact of hair cell preservation in cochlear implantation: Combined electric and acoustic hearing. *Otology & Neurotology*, 31, 1227-1232.
- Valente, C., Abrams, H. B., Benson, D., Chisolm, T. H., Citron, D., & Hampton, D. (2006). Guidelines in the audiological management of adult hearing impairment. *Audiology Today*, 18, 32-36.
- Van de Heyning, P., & Kleine Punte, A. (Eds.). (2010). *Cochlear Implants and Hearing Preservation*. (Vol. 67). Basel, Switzerland: Karger AG.
- Vermeire, K., Anderson, I., Flynn, M., Van de & Heyning, P. (2008). The influence of different speech processor and hearing aid settings on speech perception outcomes in electric acoustic stimulation patients. *Ear and Hearing*, 29(1), 76-86. doi:10.1097/AUD.1090b1013e31815d36326.
- von Ilberg, C., Kiefer, J., Tillein, J., Pfenningdorff, T., Hartmann, R., Sturzebecher, E., & Klinke, R. (1999). Electric-acoustic stimulation of the auditory system. New technology for severe hearing loss. *Journal for Oto-rhino-laryngology and Its Related Specialities*, 61, 334-340.
- Woodson, E. A., Reiss, L. A., Turner, C. W., Gfeller, K., & Gantz, B. J. (2009). The hybrid cochlear implant: A review. *Advances in Oto-rhino-laryngology*, 67, 125-134.
- Zhang, T., Dorman, M. F., & Spahr, A. J. (2010). Information from the voice fundamental frequency (F0) region accounts for the majority of the benefit when acoustic stimulation is added to electric stimulation. *Ear Hear*, 31(1), 63-69.